EUROPEAN PATENT SPECIFICATION

- (45) Date of publication and mention of the grant of the patent: 16.08.2000 Bulletin 2000/33
- (21) Application number: 93915279.9
- (22) Date of filing: 09.06.1993

- (51) Int. Cl.7: C12N 15/27, C12P 21/02, C07K 14/53, C12N 15/12, G01N 33/68
- (86) International application number: PCT/US93/05548
- (87) International publication number: WO 93/25687 (23.12.1993 Gazette 1993/30)
- (54) CRYSTALLIZATION OF M-CSF KRISTALLISIERUNG VON M-CSF CRISTALLISATION DU M-CSF
- (84) Designated Contracting States: AT BE CH DE DK ES FR GB GR IE IT LI LU MC NL PT SE
- (30) Priority: 09.06.1992 US 896512
- (43) Date of publication of application: 30.08.1995 Bulletin 1995/35
- (60) Divisional application: 99108304.9 / 0 955 365
- (73) Proprietors:
 - CHIRON CORPORATION Emeryville California 94608-2916 (US)
 - The Regents of the University of California Oakland, CA 94607-5200 (US)
- (72) Inventors:
 - PANDIT, Jayvardhan Mystic, CT 06355 (US)
 - JANCARIK, Jarmila Walnut Creek, CA 94596 (US)
 - · KIM, Sung-Hou Moraga, CA 94556 (US)
 - KOTHS, Kirston El Cerrito, CA 94530 (US)
 - · HALENBECK, Robert San Rafael, CA 94901 (US)
 - FEAR, Anna, Lisa Oakland, CA 94609 (US)
 - TAYLOR, Eric Oakland, CA 94609 (US)

- · YAMAMOTO, Ralph Martinez, CA 94553 (US)
- BOHM, Andrew Berkeley, CA 94704 (US)
- (74) Representative: Ritter, Stephen David et al Mathys & Squire 100 Grays Inn Road London WC1X 8AL (GB)
- (56) References cited: EP-A- 0 344 796

WO-A-88/08003

EP-A- 0 402 070 WO-A-92/21029

- · R. K. SCOPES 'protein purification advanced texts in chemistry' 1987, SECOND EDITION. PAGES 297-321. SPRINGER-VERLAG, NEW YORK, US.
- JOURNAL OF BIOLOGICAL CHEMISTRY. vol. 251, no. 20, 25 October 1976, BALTIMORE US pages 6300 - 6303 A. MCPHERSON, JR 'Crystallization of proteins from polyethylene glycol' cited in the application
- SCIENCE vol. 258, 20 November 1992, LANCASTER, PA US pages 1358 - 1362 J. PANDIT ET AL 'Three-dimensional structure of dimeric human recombinant Macrophage Colony-Stimulating Factor'
- · SPORN M.B. ET ROBERTS A.B. 'peptide growth factors and their receptors, handbook of experimental pharmacology vol. 1; vol. 95/1. pages 667-698.' 1990, SPINGER VERLAG, DE, BERLIN.

Note: Within nine months from the publication of the mention of the grant of the European patent, any person may give notice to the European Patent Office of opposition to the European patent granted. Notice of opposition shall be filed in a written reasoned statement. It shall not be deemed to have been filed until the opposition fee has been paid. (Art. 99(1) European Patent Convention).

- PROCEEDINGS OF THE NATIONAL ACADEMY
 OF SCIENCES OF USA. vol. 84, March 1987,
 WASHINGTON US pages 1268 1271 Y. G.
 YEUNG ET AL 'Purification of the colony
 stimulating factor 1 receptor and
 demonstration of its tyrosine kinase activity'
- JOURNAL OF MOLECULAR BIOLOGY vol. 205, no. 4, 1989, pages 783 - 785 J. M. LA LONDE ET AL 'Crystallization and preliminary X-ray studies of recombinant human granulocyte - macrophage colony stimulating factor'
- SCIENCE vol. 254, 8 November 1991, LANCASTER, PA US page 821825 B. C.
 CUNNINGHAM ET AL 'Dimerization of the extracellular domain of the human growth hormone receptor by a single hormone molecule'
- JONES, N.D. ET AL.: "Crystallization of authentic recombinant human growth hormone", BIOTECHNOLOGY, vol. 5, pages 499 to 500, May 1987

Description

10

CRYSTALLIZATION OF M-CSF

[0001] Work described herein was funded with the United States Government support. The United States Government has certain rights in inventions arising as part of that work. This application is a continuation-in-part of U.S. Serial No. 07/896,512 filed June 9, 1992.

FIELD OF THE INVENTION

[0002] The present invention relates in general to crystalline compositions of macrophage colony stimulating factor "M-CSF" and in particular to methods for the use of structural information (including X-ray diffraction patterns) of crystalline M-CSF for agonist and antagonist production, as well as assays for detection of same.

15 BACKGROUND OF THE INVENTION

[0003] Monocyte-macrophage colony-stimulating factor is produced by a variety of cells, including macrophages, endothelial cells and fibroblasts (see, Ralph et al., "The Molecular and Biological Properties of the Human and Murine Members of the CSF-1 Family" in Molecular Basis of Lymphokine Action, Humana Press, Inc., (1987),). M-CSF is composed of two "monomer" polypeptides, which form a biologically active dimeric M-CSF protein (hereinafter referred to as "M-CSF dimer"). M-CSF belongs to a group of biological agonists that promote the production of blood cells. Specifically, it acts as a growth and differentiation factor for bone marrow progenitor cells of the mononuclear phagocyte lineage. Further, M-CSF stimulates the proliferation and function of mature macrophages via specific receptors on responding cells. In clinical trials M-CSF has shown promise as a pharmaceutical agent in the correction of blood cell deficiencies arising as a side-effect of chemotherapy or radiation therapy for cancer and may be beneficial in treating fungal infections associated with bone marrow transplants. M-CSF may also play significant biological roles in pregnancy, uveitis, and atherosclerosis. Development of M-CSF agonists or antagonists may prove to be of value in modifying the biological events involved in these conditions.

[0004] M-CSF exists in at least three mature forms: short (M-CSFα), intermediate (M-CSFγ), and long (M-CSFβ). Mature M-CSF is defined as including polypeptide sequences contained within secreted M-CSF following amino terminus processing to remove leader sequences and carboxyl terminus processing to remove domains including a putative transmembrane region. The variations in the three mature forms are due to alternative mRNA splicing (see, Cerretti et al. Molecular Immunology, 25:761 (1988)). The three forms of M-CSF are translated from different mRNA precursors, which encode polypeptide monomers of 256 to 554 amino acids, having a 32 amino acid signal sequence at the amino terminal and a putative transmembrane region of approximately 23 amino acids near the carboxyl terminal. The precursor peptides are subsequently processed by amino terminal and carboxyl terminal proteolytic cleavages to release mature M-CSF. Residues 1-149 of all three mature forms of M-CSF are identical and are believed to contain sequences essential for biological activity of M-CSF. In vivo M-CSF monomers are dimerized via disulfide-linkage and are glycosylated. Some, if not all, forms of M-CSF can be recovered in membrane-associated form. Such membrane-bound M-CSF may be cleaved to release secreted M-CSF. Membrane associated M-CSF is believed to have biological activity similar to M-CSF, but may have other activities including cell-cell association or activation.

[0005] Polypeptides, including the M-CSFs, have a three-dimensional structure determined by the primary amino acid sequence and the environment surrounding the polypeptide. This three-dimensional structure establishes the polypeptide's activity, stability, binding affinity, binding specificity, and other biochemical attributes. Thus, a knowledge of a protein's three-dimensional structure can provide much guidance in designing agents that mimic, inhibit, or improve its biological activity in soluble or membrane bound forms.

[0006] The three-dimensional structure of a polypeptide may be determined in a number of ways. Many of the most precise methods employ X-ray crystallography (for a general review, see, Van Holde, *Physical Biochemistry, Prentice-Hall, N.J.* pp. 221-239, (1971), which is incorporated herein by reference). This technique relies on the ability of crystalline lattices to diffract X-rays or other forms of radiation. Diffraction experiments suitable for determining the three-dimensional structure of macromolecules typically require high-quality crystals. Unfortunately, such crystals have been unavailable for M-CSF as well as many other proteins of interest. Thus, high-quality, diffracting crystals of M-CSF would assist the determination of its three-dimensional structure.

[0007] Various methods for preparing crystalline proteins and polypeptides are known in the art (see, for example, McPherson, et al. "Preparation and Analysis of Protein Crystals", A. McPherson, Robert E. Krieger Publishing Company, Malabar, Florida (1989); Weber, Advances in Protein Chemistry 41:1-36 (1991); U.S. Patent No. 4,672,108; and U.S. Patent No. 4,833,233;). Although there are multiple approaches to crystallizing polypeptides, no single set of conditions provides a reasonable expectation of success, especially when the crystals must be suitable for X-ray dif-

fraction studies. Thus, in spite of significant research, many proteins remain uncrystallized.

[0008] In addition to providing structural information, crystalline polypeptides provide other advantages. For example, the crystallization process itself further purifies the polypeptide, and satisfies one of the classical criteria for homogeneity. In fact, crystallization frequently provides unparalleled purification quality, removing impurities that are not removed by other purification methods such as HPLC, dialysis, conventional column chromatography, etc. Moreover, crystalline polypeptides are often stable at ambient temperatures and free of protease contamination and other degradation associated with solution storage. Crystalline polypeptides may also be useful as pharmaceutical preparations. Finally, crystallization techniques in general are largely free of problems such as denaturation associated with other stabilization methods (e.g. lyophilization). Thus, there exists a significant need for preparing M-CSF compositions in crystalline form and determining their three-dimensional structure. The present invention fulfills this and other needs. Once crystallization has been accomplished, crystallographic data provides useful structural information which may assist the design of peptides that may serve as agonists or antagonists. In addition, the crystal structure provides information useful to map, the receptor binding domain which could then be mimicked by a small non-peptide molecule which may serve as an antagonist or agonist.

SUMMARY OF THE INVENTION

15

[0009] The present invention provides crystalline forms of M-CSF dimers. Preferably, the dimers are formed from polypeptides containing between 146 to 162 amino acids residues at or near the N-terminus of mature M-CSF (e.g. glu_1 glu_2 $val_3...$). In a specific embodiment, the polypeptide includes residues 4 to 158 of mature M-CSF α polypeptide, preferably in the non-glycosylated form.

[0010] Another aspect of the invention provides a method of crystallizing an M-CSF. A preferred crystallization method according to the present invention includes the following steps: mixing a preselected, substantially pure M-CSF dimer and a precipitant to form an M-CSF mixture; precipitating crystalline M-CSF from the mixture; and isolating the M-CSF in crystalline form. In some specific embodiments, the precipitant contains polyethylene glycol. Other components such as ammonium sulfate and/or other ionic compounds may be added to the solution. It has been found by x-ray crystallography that M-CSF produced by the method of the present invention can crystallize into the P2₁2₁2₁ space groups for example.

[0011] Variations of the crystallization method are also provided. For example, the step of precipitating crystals from the M-CSF mixture may involve equilibrating the M-CSF mixture with a second mixture. The second mixture is typically a solution that consists of a higher concentration of precipitant than the first M-CSF mixture. The step of equilibrating preferably consists of applying the M-CSF mixture to a surface and allowing the applied M-CSF mixture to come into equilibrium with a reservoir of the second mixture. In other embodiments, the step of precipitating M-CSF crystals is initiated by seeding the M-CSF mixture with seed crystals or altering the temperature of the M-CSF mixture. Another aspect of the invention involves identifying compounds that have structures that mimic a receptor binding region of the three-dimensional structure of M-CSF to varying degrees and can in many instances function as M-CSF agonists or antagonists. Compounds that interact with the receptor-binding region of M-CSF may be antagonists. The three-dimensional alpha-carbon coordinates of a truncated M-CSF dimer is presented in Appendix 1. In one embodiment of the present invention, the three-dimensional structure of M-CSF is obtained by first crystallizing an M-CSF dimer (having M-CSF receptor-binding residues) to form at least one M-CSF crystal. Next, a source of radiation is used for irradiating an M-CSF crystal to obtain a diffraction pattern of the M-CSF crystal. Finally, a three-dimensional structure of M-CSF is obtained from the diffraction pattern. In most embodiments, the three-dimensional structure includes an M-CSF receptor-binding region.

Structural features identified in accordance with the invention may be used for selecting candidate amino acid substitutions in a protein, based on structural information, and more particularly M-CSF, comprising determining the three-dimensional structure of M-CSF by the methods of the present invention; followed by determining the solvent accessible amino acid residues of the protein, determining which residues are not involved in dimer formation. Applying these criteria, amino acids in M-CSF which are solvent accessible and which are not involved in dimer formation are selected for substitution with non-conservative amino acids. Since M-CSF β has intrachain disulfide bonds involving cysteines 157 and/or 159, we believe the C-terminal region of M-CSF extends from the "rear" of the structure we have solved, providing a variable-length "tether" for membrane-bound forms of M-CSF. Thus, the "front" or receptor-binding region of M-CSF is on the opposite side of the molecules, consisting of solvent-accessible residues in or near helices A, C, and D, including residues from about 6 to 26, 71 to 90, and 110 to 130, respectively, of native M-CSF. Preferred amino acids for substitution and preferred substituting amino acids include but are not limited to: H15 \rightarrow A or L; Q79 \rightarrow A or D; R86 \rightarrow E or D; E115 \rightarrow A; E41 \rightarrow K or R; K93 \rightarrow A or E; D99 \rightarrow K or R; L55 \rightarrow Q or N; S18 \rightarrow A or K; Q20 \rightarrow A or D; 175 \rightarrow K or E; V78 \rightarrow K or R; L85 \rightarrow E or N; D69 \rightarrow K or R; N70 \rightarrow A or E; H9 \rightarrow A or D; N63 \rightarrow K or R; and T34 \rightarrow Q or K. Most preferred are those substitutions give rise to novel M-CSF agonists and M-CSF antagonists. Additionally, the present invention is also directed to a method for producing antagonists and agonists by substituting at least

one and preferably fewer than 5 solvent accessible residues per M-CSF monomer.

[0013] The invention also enables the production of heterodimeric M-CSF in which only one subunit contains substituted solvent accessible amino acids involved in signal transduction and to heterodimeric M-CSF in which each subunit contains different substituted solvent accessible amino acids involved in signal transduction. The present invention is also enables the production of M-CSF having amino acid substitutions which do not impair binding to the M-CSF receptor. Screening for agonists and antagonists is then accomplished using bioassays and receptor binding assays using methods well known in the art, including those described in the Examples below.

[0014] In addition the invention is directed to an isolated, purified, soluble and functional M-CSF receptor. The present invention is also directed to a method for screening M-CSF agonists and antagonists using a soluble M-CSF receptor.

[0015] A further understanding of the present invention can be obtained by reference to the drawings and discussion of specific embodiments.

BRIEF DESCRIPTION OF THE FIGURES

[0016]

15

20

30

35

40

Figure 1 is a section of a diffraction pattern of an M-CSFα crystal prepared according to the present invention; Figure 2 is a topology diagram showing the disulfide bonds in truncated dimeric M-CSF;

Figure 3 is a stereodiagram of the C-alpha backbone with every tenth residue labelled and with the non-crystallographic symmetry axis indicated by a dotted line;

Figures 4A and 4B present two views of a ribbon diagram highlighting the secondary structural elements of M-CSF. The cysteine residues have been represented by a ball-and stick model and the non-crystallographic symmetry axis is indicated by a dotted line;

Figures 5A-C illustrate size-exclusion HPLC analysis of the NΔ3CΔ158 M-CSF short clone homodimer (158) (Figure 5A), NΔ3CΔ221 C157S, C159S long clone homodimer (221F) (Figure 5B), and the short clone/long clone heterodimer (158/221F) (Figure 5C) and their corresponding biological activities;

Figures 6A-6C illustrates size exclusion HPLC (Figure 6A) and both non-reduced (Figure 6B) and reduced (Figure 6C) SDS-PAGE analysis of the preparative purification of M-CSF, more particularly, Figure 6A graphically illustrates the separation on Phenyl-HPLC size exclusion chromatography of the three species of M-CSF dimers of Figures 5A-5C, i.e., the 158 homodimer, the 221F homodimer and the 158/221 heterodimer, and indicates that the three absorbance peaks at 280 nm (solid line) correlate with M-CSF activity in U/ml x 10⁻⁶ (dotted line); Figure 6B illustrates an SDS-PAGE analysis under non-reducing conditions of preparative purification of the 158/221F heterodimer (intermediate molecular weight species) relative to the 158 homodimer (lower molecular weight species) and the 221F homodimer (highest molecular weight species); Figure 6C illustrates an SDS-PAGE analysis under reducing conditions of the preparative purification of the 158/221F heterodimer (middle lanes) relative to the 158 homodimer (left lanes) and the 221F homodimer (right lanes); and

Figure 7 illustrates the competitive binding of M-CSF and M-CSF muteins to NFS60 cell M-CSF receptors. In Figure 7, competitive binding curves are shown for M-CSFα NΔ3CΔ158 (closed circles); M-CSFα NΔ3CΔ158 H9A, H15AM-CSFβ NΔ3CΔ221 C157S, C159S heterodimer (closed squares); dimeric Q20A, V78KF mutein (open circles); and dimeric H9A, H15A mutein (open squares).

DESCRIPTION OF THE PREFERRED EMBODIMENTS DEFINITIONS

45 [0017] As used herein "M-CSF polypeptide" refers to a human polypeptide having substantially the same amino acid sequence as the mature human M-CSFα, M-CSFβ, or M-CSFβ polypeptides described in Kawasaki et al, Science 230:291 (1985), Cerretti et al, Molecular Immunology, 25:761 (1988), or Ladner et al., EMBO Journal 6:2693 (1987). Such terminology reflects the understanding that the three mature M-CSFs have different amino acid sequences, as described above.

[0018] Certain modifications to the primary sequence of M-CSF can be made by deletion, addition, or alteration of the amino acids encoded by the DNA sequence without destroying the desired structure (e.g., the receptor binding ability of M-CSF) in accordance with well-known recombinant DNA techniques. Further, a skilled artisan will appreciate that individual amino acids may be substituted or modified by oxidation, reduction or other derivitization, and the polypeptide may be cleaved to obtain fragments that retain the active binding site and structural information. Such substitutions and alterations result in polypeptides having an amino acid sequence which falls within the definition of polypeptide "having substantially the same amino acid sequence as the mature M-CSFα, M-CSFβ, and M-CSFγ polypeptides."

[0019] For purposes of crystallization, preferred lengths of the M-CSF α , β or γ monomers are between about 145 and 180 amino acids (counting from the mature amino terminus), and more preferably between about 145 and 162

amino acids long. A specific monomer that may be present in a crystallizable dimer is M-CSF α and is N Δ 3M-CSF α C Δ 158 (3 amino acids are deleted from the amino terminus and the total length is 155 amino acids). All lengths are inclusive. As used herein the term "M-CSF α (4-158)" denotes an M-CSF having amino acid residues 4 to 158 of the mature, processed M-CSF α polypeptide. Other nomenclature designations for C-terminal and N-terminal truncations of native M-CSF are set forth in U.S. Patent No. 4,929,700.

[0020] Crystallizable glycosylation variants of the M-CSF polypeptides are included within the scope of this invention. These variants include polypeptides completely lacking in glycosylation and variants having at least one fewer glycosylated site than the mature forms, as well as variants in which the glycosylation pattern has been changed from the native forms. Also included are deglycosylated and unglycosylated amino acid sequence variants, as well as deglycosylated and unglycosylated M-CSF subunit having the mature amino acid sequence (see, U.S. Patent No. 5,032,626). [0021] "M-CSF" dimer refers to two M-CSF polypeptide monomers that have dimerized M-CSF dimers may include two identical polypeptide monomers (homodimers) or two different polypeptide monomers (heterodimers such as an M-CSFα-M-CSFβ dimer, an M-CSF long clone and short clone dimer). M-CSF monomers may be convened to M-CSF dimers *in vitro* as described in U.S. Patent No. 4,929,700. Recombinantly expressed M-CSFs may also be variably glycosylated as they exist *in vivo*, partially glycosylated, or completely lacking in glycosylation (unglycosylated). Glycosylated M-CSFs may be produced *in vivo* with carbohydrate chains which may later be enzymatically deglycosylated *in vitro*.

[0022] Biologically active M-CSF exhibits a spectrum of activity understood in the art. For instance, M-CSF stimulates the proliferation and function of mature macrophages via specific receptors on responding cells. Further, M-CSF acts as a mononuclear phagocyte progemtor growth factor. The standard *in vitro* colony stimulating assay of Metcalf, *J. Cell Physiol.* 76:89 (1970) results primarily in the formation of macrophage colonies when M-CSF is applied to stem cells. Other biological assays are based on M-CSF induced proliferation of M-CSF dependent cells such as the NFS-60 cell line. As used herein "M-CSF having biological activity" refers to M-CSF, including fragments and sequence variants thereof as described herein; that exhibit an an-recognized spectrum of activity with respect to biological systems. Such M-CSF having biological activity will typically have certain structural attributes in common with those of the mature M-CSF forms such as receptor binding site tertiary structure.

[0023] Agonists are substances that exhibit greater activity *per se* than the native ligand while antagonists are substances that suppress, inhibit, or interfere with the biological activity of a native ligand. Agonists and antagonists may be produced by the methods of the present invention for use in the treatment of diseases in which M-CSF has been implicated either as a potential treatment (e.g., for treating blood cell deficiencies arising as a side effect of chemotherapy treating fungal infection associated with bone marrow transplants and others) or as having a role in the pathogenesis of the disease (e.g., ovarian cancer, uveitis, atherosclerosis).

[0024] Crystallization of M-CSF species in accordance with the present invention includes four general steps: expression, purification, crystallization and isolation.

Expression of Recombinant M-CSF

35

[0025] M-CSF crystallization requires an abundant source of M-CSF that may be isolated in a relatively homogeneous form. A variety of expression systems and hosts are suitable for the expression of M-CSF and will be readily apparent to one of skill in the art. Because of the variability of glycosylation and other post-transnational modifications present in M-CSF produced in certain eukaryotic hosts, expression in *E. coli* may provide M-CSF with advantageous properties with regard to crystallization. Typical *in vitro* M-CSF expression systems are described in U.S. Patent No. 4,929,700, for example.

[0026] For use in the present invention, a variety of M-CSF polypeptides can also be readily designed and manufactured utilizing recombinant DNA techniques well known to those skilled in the art. For example, the M-CSF amino acid sequence can vary from the naturally occurring sequence at the primary structure level by amino acid substitutions, insertions, deletions, and the like. These modifications can be used in a number of combinations to produce the final modified polypeptide chain. The present invention is useful for crystallizing such polypeptides and dimers thereof. [0027] In general, modifications of the genes encoding the M-CSF polypeptide are readily accomplished by a variety of well-known techniques, such as site-directed mutagenesis (see, Gillman and Smith, Gene 8:81-97 (1979) and Roberts, S. et al., Nature 328:731-734 (1987) and U.S. Patent No. 5,032,676.). Most modifications are evaluated by screening in a suitable assay for the desired characteristic. For instance, a change in the M-CSF receptor-binding character of the polypeptide can be detected by competitive assays with an appropriate reference polypeptides or by the bioassays described in U.S. Patent No. 4,847,201, which is incorporated herein by reference.

[0028] Insertional variants of the present invention are those in which one or more amino acid residues are introduced into a predetermined site in the M-CSF. For instance, insertional variants can be fusions of heterologous proteins or polypeptides to the amino or carboxyl terminus of the subunits. Substitutional variants are those in which at least one residue has been removed and a different residue inserted in its place. Nonnatural amino acids (i.e., amino acids not

normally found in native proteins), as well as isosteric analogs (amino acid or otherwise) are also suitable for use in this invention. Examples of suitable substitutions are well known in the art, such as the Glu->Asp, Ser->Cys, and Cys->Ser, His-> alanine for example. Another class of variants are deletional variants, which are characterized by the removal of one or more amino acid residues from the M-CSF.

[0029] Other variants of the present invention may be produced by chemically modifying amino acids of the native protein (e.g., diethylpyrocarbonate treatment which modifies histidine residues). Preferred or chemical modifications which are specific for certain amino acid side chains. Specificity may also be achieved by blocking other side chains with antibodies directed to the side chains to be protected. Chemical modification includes such reactions as oxidation, reduction, amidation, deamidation, or substitution of bulky groups such as polysaccharides or polyethylene glycol (see e.g., U.S. Patent No. 4,179,337 and WO91/21029).

[0030] Exemplary modifications include the modification of lysinyl and amino terminal residues by reaction with succinic or other carboxylic acid anhydrides. Modification with these agents has the effect of reversing the charge of the lysinyl residues. Other suitable reagents for modifying amino-containing residues include imidoesters such as methyl picolinimidate; pyridoxal phosphate; pyridoxal chloroborohydride; trinitrobenzenesulfonic acid; O-methylisourea, 2,4-pentanedione; and transaminase-catalyzed reaction with glyoxylate, and N-hydroxysuccinamide esters of polyethylenene glycol or other bulky substitutions.

[0031] Arginyl residues may be modified by reaction with a number of reagents, including phenylglyoxal, 2,3-butanedione, 1,2-cyclohexanedione, and ninhydrin. Modification of arginine residues requires that the reaction be performed in alkaline conditions because of the high pK_a of the guanidine functional group. Furthermore, these reagents may react with the groups of lysine as well as the arginine epsilon-amino group.

[0032] Tyrosyl residues may also be modified with particular interest in introducing spectral labels into tyrosyl residues by reaction with aromatic diazonium compounds or tetranitromethane are used to form O-acetyl tyrosyl species and 3-nitro derivatives, respectively. Tyrosyl residues may also be iodinated using ¹²⁵I or ¹³¹I to prepare labeled proteins for use in radioimmunoassay.

[0033] Carboxyl side groups (aspartyl or glutamyl) may be selectively modified by reaction with carbodiimides (R-N=C=N-R¹), where R and R¹ are different alkyl groups, such as 1-cyclohexyl-3-(2-morpholinyl-4-ethyl) carbodiimide or 1-ethyl-3-(4-azonia-4,4-dimethylpentyl) carbodiimide. Furthermore, aspartyl and glutamyl residues are convened to asparaginyl and glutaminyl residues by reaction with ammonium ions.

[0034] Conversely, glutaminyl and asparaginyl residues may be deamidated to the corresponding glutamyl and aspartyl residues, respectively, under mildly acidic conditions. Either form of these residues falls within the scope of this invention.

[0035] Other modifications include hydroxylation of proline and lysine, phosphorylation of hydroxyl groups of seryl or threonyl residues, methylation of the α -amino groups of lysine, arginine, and histidine side chains (T.E. Creighton, *Proteins: Structure and Molecular Properties*, W.H. Freeman & Co., San Francisco, pp. 79-86 [1983]), acetylation of the N-terminal amine, and amidation of any C-terminal carboxyl group.

[0036] The availability of a DNA sequence encoding M-CSF permits the use of various expression systems to produce the desired polypeptides. Construction of expression vectors and recombinant production from the appropriate DNA sequences are performed by methods well known in the art. These techniques and various other techniques are generally performed according to Sambrook et al., Molecular Cloning - A Laboratory Manual, Cold Spring Harbor Laboratory, Cold Spring Harbor, New York (1989), and Kriegler, M., Gene Transfer and Expression, A Laboratory Manual, Stockton Press, New York (1990).

Purification of M-CSF

10037] Purification steps are employed to ensure that the M-CSF is isolated, prior to crystallization, in a relatively homogeneous state. In general, a higher purity solution increases the likelihood of success of subsequent crystallization steps. Typical purification methods include the use of centrifugation, partial fractionation using salt or organic compounds, dialysis, conventional column chromatography (such as ion-exchange, molecular sizing chromatography etc.), high performance liquid chromatography (HPLC), and gel electrophoresis methods (see, e.g., Deutcher, "Guide to Protein Purification" in Methods in Enzymology (1990), Academic Press, Berkely, CA, which is incorporated herein by reference for all purposes). Preferred purification conditions for generating unusually homogeneous M-CSF species as well as purification of these species are disclosed, for example, in U.S. Patent No. 4,929 700. Other purification methods are known and will be apparent to one of skill in the art.

Crystallization of M-CSF

[0038] Although many of the same physical principles govern crystallization of polypeptides (including M-CSF dimers) and small molecules, the actual crystallization mechanisms differ significantly. For example, the lattice of small-

molecule crystals effectively excludes solvent while that of polypeptide crystals includes substantial numbers of solvent molecules. Thus, special techniques must typically be applied to crystallize polypeptides.

[0039] Polypeptide crystallization occurs in solutions where the polypeptide concentration exceeds its solubility maximum (i.e., the polypeptide solution is supersaturated). Such "thermodynamically metastable" solutions may be restored to equilibrium by reducing the polypeptide concentration, preferably through precipitation of the polypeptide crystals. Often polypeptides may be induced to crystallize from supersaturated solutions by adding agents that alter the polypeptide surface charges or perturb the interactions between the polypeptide and bulk water to promote associations that lead to crystallization.

[0040] Compounds known as "precipitants" are often used to decrease the solubility of the polypeptide in a concentrated solution. Precipitants induce crystallization by forming an energetically unfavorable precipitant-depleted layer around the polypeptide molecules. To minimize the relative amount of this depletion layer, the polypeptides form associations and ultimately crystals as explained in Weber, Advances in Protein Chemistry 41:1-36 (1991). In addition to precipitants, other materials are sometimes added to the polypeptide crystallization solution. These include buffers to adjust the pH of the solution (and hence surface charge on the peptide) and salts to reduce the solubility of the polypeptide. Various precipitants are known in the art and include the following: ethanol, 3-ethyl-2,4 pentanediol; and many of the polyglycols, such as polyethylene glycol. A suitable precipitant for crystallizing M-CSF is polyethylene glycol (PEG), which combines some of the characteristics of the salts and other organic precipitants (see, for example, Ward et al., J. Mol. Biol. 98:161 [1975] and McPherson J. Biol. Chem. 251:6300 [1976],).

[0041] Commonly used polypeptide crystallization methods include the following techniques: batch, hanging drop, seed initiation, and dialysis. In each of these methods, it is important to promote continued crystallization after nucleation by maintaining a supersaturated solution. In the batch method, polypeptide is mixed with precipitants to achieve supersaturation, the vessel is sealed and set aside until crystals appear. In the dialysis method, polypeptide is retained in a sealed dialysis membrane which is placed into a solution containing precipitant. Equilibration across the membrane increases the polypeptide and precipitant concentrations thereby causing the polypeptide to reach supersaturation levels.

25

[0042] In the hanging drop technique, an initial polypeptide mixture is created by adding a precipitant to concentrated polypeptide solution. The concentrations of the polypeptide and precipitants are such that in this initial form, the polypeptide does not crystallize. A small drop of this mixture is placed on a glass slide which is inverted and suspended over a reservoir of a second solution. The system is then sealed. Typically the second solution contains a higher concentration of precipitant or other dehydrating agent. The difference in the precipitant concentrations causes the protein solution to have a higher vapor pressure than the solution. Since the system containing the two solutions is sealed, an equilibrium is established, and water from the polypeptide mixture transfers to the second solution. This equilibration increases the polypeptide and precipitant concentration in the polypeptide solution. At the critical concentration of polypeptide and precipitant, a crystal of the polypeptide will form. The hanging drop method is well known in the art (see, McPherson J. Biol. Chem. 251:6300 [1976],).

[0043] Another method of crystallization introduces a nucleation site into a concentrated polypeptide solution. Generally, a concentrated polypeptide solution is prepared and a seed crystal of the polypeptide is introduced into this solution. If the concentrations of the polypeptide and any precipitants are correct, the seed crystal will provide a nucleation site around which a larger crystal forms.

In preferred embodiments, the crystals of the present invention will be formed from a dimer of M-CSF polypeptides. Preferred crystals are typically at least about 0.2 x 0.2 x 0.05 mm, more preferably larger than 0.4 x 0.4 x 0.4 mm, and most preferably larger than 0.5 x 0.5 x 0.5 mm. After crystallization, the protein may be separated from the crystallization mixture by standard techniques.

[0045] The crystals so produced have a wide range of uses. For example, high quality crystals are suitable for X-ray or neutron diffraction analysis to determine the three-dimensional structure of the M-CSF and, in particular, to assist in the identification of its receptor binding site. Knowledge of the binding site region and solvent-accessible residues available for contact with the M-CSF receptor allows rational design and construction of agonists and antagonist for M-CSFs. Crystallization and structural determination of M-CSF muteins having altered receptor binding ability or bioactivity allows the evaluation of whether such changes are caused by general structural deformation or by side chain alteration at the substitution site.

[0046] In addition, crystallization itself can be used as purification method. In some instances, a polypeptide or protein will crystallize from a heterogeneous mixture into crystals. Isolation of such crystals by filtration, centrifugation, etc. followed by redissolving the polypeptide affords a purified solution suitable for use in growing the high-quality crystals necessary for diffraction studies. These high-quality crystals may also be dissolved in water and then formulated to provide an aqueous M-CSF solution having various uses known in the art including pharmaceutical purposes.

[0047] Of course, amino acid sequence variants of M-CSF may also be crystallized and used. These mutants can be used for, among other purposes, obtaining structural information useful for directing modification of the binding affinity for M-CSF receptors. As with the naturally occurring forms, the modified M-CSF forms may be useful as pharmaceu-

tical agents for stimulating bone marrow proliferation, overcoming immune suppression and fungal diseases induced by chemotherapy, improving therapeutic efficacy, and lessening the severity or occurrence of side effects during therapeutic use of the present invention. Furthermore, modified M-CSFs may be useful for treatment of disease in which soluble or membrane-bound M-CSF causes or exacerbates the disease state.

Characterization of M-CSF

5

[0048] After purification, crystallization and isolation, the subject crystals may be analyzed by techniques known in the art. Typical analysis yield structural, physical, and mechanistic information about the peptides. As discussed above, X-ray crystallography provides detailed structural information which may be used in conjunction with widely available molecular modeling programs to arrive at the three-dimensional arrangement of atoms in the peptide. Exemplary modeling programs include "Homology" by Biosym (San Diego, California), "Biograf" by BioDesign, "Nemesis" by Oxford Molecular, "SYBYL" and "Composer" by Tripos Associates, "CHARM" by Polygen (Waltham, Massachusetts), "AMBER" by University of California, San Francisco, and "MM2" and "MMP2" by Molecular Design, Ltd.

[0049] Peptide modeling can be used to design a variety of agents capable of modifying the activity of the subject peptide. For example, using the three-dimensional structure of the active site, agonists and antagonists having complementary structures can be designed to enhance the therapeutic utility of M-CSF treatment or to block the biological activity of M-CSF. Further, M-CSF structural information is useful for directing design of proteinaceous or non-proteinaceous M-CSF agonists and antagonists, based on knowledge of the contact residues between the M-CSF ligand and its receptor. Such residues are identified by the M-CSF crystal structure as those which are solvent-accessible, distal to the carboxyl terminal membrane anchoring region not involved in dimer interface stabilizations, and possibly including residues not conserved between human and mouse M-CSF (which does not recognize the human M-CSF receptor).

25 Example 1

[0050] Systematic crystallization trials with M-CSF were made using the hanging drop technique. A microdroplet (5µl) of mother liquor is suspended from the underside of a microscope cover slip, which is placed over a well containing 1ml of the precipitating solution. 60-70 initial trials were set up, in which pH, temperature, counterion and precipitant were varied. From these trials, the few that gave promising microcrystals were picked for more careful examination.

[0051] It was discovered that suitable crystals may be grown from a 20 μ l drop containing: 10 mg/ml protein, 100 mM MgCl₂, 50 mM Tris.Cl, pH 8.5, and 12% PEG 4000. This drop was equilibrated against a reservoir containing 24% PEG 4000. Tiny, needle-like crystals appeared in 2-3 days which were redissolved in 10 μ l water and recrystallized at room temperature. Good quality chunky crystals appeared in 7-9 days in sizes ranging from 0.3 x 0.3 x 0.3 mm to 0.5 x 0.5 x 1.0 mm.

[0052] Precession photographs revealed the space group to be P2₁2₁2₁ with unit cell dimensions: a=33.53Å, b=65.11Å, c=159.77Å. This gives a unit cell volume of 349084.5Å³, which is consistent with a dimer in the crystallographic asymmetric unit, and 52% of the unit cell volume being occupied by solvent. The crystals diffracted to a resolution of 3Å on a Rigaku rotating anode X-ray generator (Danvers, Massachusetts) operated at 50 kV and 60mA, and to 2.6Å in synchrotron radiation.

[0053] Screening for heavy atom derivatives was done by soaking crystals into solutions of heavy-metal salts. Zero-level precession pictures of the soaks were used to identify potential derivatives. About 30 different soaking conditions were examined, of which 4 potential derivatives were identified. Unfortunately, some soaks caused the crystals to exhibit non-isomorphism (i.e., the heavy atom soaks induced a change in cell dimensions, making them unusable for phase calculation).

[0054] Three-dimensional intensity data were collected on film using an oscillation camera on the X-ray beam-line at the National Synchrotron Light Source, Brookhaven. Several other data sets, of native (underivatized M-CSF) as well as potential derivative crystals have been collected on a Rigaku X-ray generator. The following data sets were collected.

Crystal	Resolution(Å)	N(observations)	N(unique)	X-Ray Source
Native	2.8	27922	7311	Synchrontron
Native	2.9	35236	7002	Rigaku (film)
Native	3.5	5144	5116	Rigaku (diffractometer)

55

50

(continued)

Crystal	Resolution(A)	N(observations)	N(unique)	X-Ray Source
K ₂ Hg(SCN) ₄	3.5	15885	4119	Rigaku (film)
UO ₂ Cl ₂	3.5	25492	5048	Rigaku (film)
Cis-Pd(NH ₃) ₂ Cl ₂	3.1	26122	6304	Synchrotron

Example 2

10

[0055] Recombinant M-CSF polypeptides were purified from $E.\ coli$. and renatured to form a disulfide-linked dimeric protein as described in U.S. Patent No. 4,929,700. Crystallization of the resulting unglycosylated M-CSF α protein (amino acids 4-158 in homodimeric form) was performed by the hanging drop method. Glass microscope plates were siliconized prior to use by dipping immersion into a 1% (volume:volume) solution of the organosilane compound, Prosil-28 (PCR Incorporated, Gainesville, Florida, 32602) washing the treated glass plates with water, and baking at 180 degrees.

[0056] A 2 mg/ml aqueous solution of purified human recombinant M-CSF was dialyzed and concentrated against 50 mM Tris-HCl (pH 8.5) using a dialysis tubing having a 10kD cutoff. The final concentration of polypeptide (10 mg/ml) was determined by ultraviolet spectrophotometry at 280 nm.

20 [0057] About 7 microliters of the concentrated solution was mixed in each well of the spot plate with 7 microliter of 20% (v/v) PEG 4000, 0.2 M MgCl₂, 0.1 M Tris-HCl (pH 8.5). The spot plate was then placed in a clear plastic sandwich box containing 20 ml of 23% PEG 4000, 0.2 M MgCl₂, 0.1 M Tris-HCl (pH 8.5) and the box was immediately sealed and stored at room temperature. Minor variations in this procedure such as altering buffer conditions are within the scope of the present invention. For example, in a preferred embodiment of the present invention, buffer conditions were altered to include 150 mM MgCl₂ and 24% PEG 4000.

[0058] After 3-5 days, small microcrystals having a size of 0.1 x 0.1 x 0.05 mm appeared in each well. These microcrystals were isolated and redissolved in 25 microliter of 50 mM Tris-HCl and allowed to stand at room temperature. The purified M-CSF crystallized from solution into large hexagonal prism shaped crystals ranging in size from 0.3 x 0.3 x 0.3 mm to 1 mm x 2 mm x 0.5 mm. These crystals were stable at room temperature for at least three months. In some instances, an artificial mother liquor was prepared using 23% PEG 4000 and 150 mM MgCl₂ crystals were then added to this mother liquor. In these cases, the crystals were removed from the mother liquor immediately prior to analysis.

[0059] Using reducing and non-reducing SDS-PAGE analysis, the M-CSF in the crystals was shown to be identical

[0059] Using reducing and non-reducing SDS-PAGE analysis, the M-CSF in the crystals was shown to be identical in molecular weight to the biologically active starting material. Thus, the M-CSF structure obtained from the crystals is likely to be essentially identical to the structure of biologically active M-CSF.

Example 3

[0060] Glass microscope slides were prepared as described in Example 2. 7 microliters of the same concentrated M-CSFα protein solution was mixed in each well of the spot plate with 7 microliter of 30% (v/v) PEG 4000, 0.2 M ammonium acetate, 0.1 M acetate buffer (pH 7.5). The spot plate was then placed in a clear plastic sandwich box containing 20 ml of 30% PEG 4000, 0.2 M ammonium acetate, 0.1 M acetate buffer (pH 7.5) and the box was immediately sealed and stored at room temperature. After 3-5 days, thin, plate-like, fragile crystals having a size of approximately 0.3 x 0.3 x 0.05 mm appeared.

45 Example 4

Preliminary X-ray Analysis

[0061] X-ray crystallographic analysis using precession photographs showed that the crystals produced in Example 2 have an orthorhombic crystal lattice in the P2₁2₁2₁ space group with cell dimensions a=33.54, b=65.26, c=159.63 d = 90.0, c = 90.0 and f = 90.0 angstroms and diffract to a nominal resolution of 2.6 angstroms using synchrotron radiation. These data provided a unit cell volume of 348084.5 angstroms³, which is consistent with a dimer in the crystallographic asymmetric unit with 52% of the unit cell volume being occupied by solvent. Fig. 1 is a 12-degree precession photograph of the Okl-zone section of the M-CSF crystal. The photograph was taken using a precession camera manufactured by Enraf-Nonius Company (Delft, Holland), mounted on a Rigaku RU-200 X-ray generator operated at 50kV and 50mA.

Example 5

Testing of M-CSF Receptor Binding Ability Using Soluble M-CSF Receptor

An essential step in the biological function of M-CSF in vivo is the binding to the M-CSF receptor, also referred to as the c-fms gene product. Recombinant human soluble M-CSF receptor (rhsM-CSFR), representing amino acids 20 to 511 (Coussens, L et al, Nature, 320:277 (1986)) was used as an in vitro assay reagent to test the receptorbinding ability of M-CSF proteins. To generate a soluble form of the transmembrane receptor, only the extracellular domain of the human M-CSF receptor was expressed in a baculovirus/insect cell recombinant expression system. In order to purify the soluble receptor without adversely affecting tertiary or quaternary structure, non-denaturing chromatographic methods were chosen, as described below. Other choices exist for the purification of the recombinant receptor. Affinity chromatography may be employed when either a suitable antibody to or ligand for the receptor are available. Alternatively, "tags" may be added to the C-terminus of the recombinant receptor, i.e., KT3 antibody recognition sequence, and purified by an anti-tag antibody, i.e., KT3, column, for use in affinity chromatography. In expression systems in which the rhsM-CSFR is glycosylated, lectin chromatography can be used to enrich for specific glycoproteins. The rhsM-CSFR can be used to study ligand/receptor interactions as well as ligand-induced receptor dimerization. The assay used to detect ligand/receptor binding employed the use of size exclusion-HPLC, essentially as described in European Patent Application WO92/21029, C. Cunningham, et al., with the following modifications: the column used was a Superose 6 (Pharmacia LKB Biotechnology, Inc.) and the mobile phase was PBS at 0.5 ml/min and a M-CSF to rhsM-CSFR ratio of 1:2. At this ratio, the M-CSF/rhsM-CSFR complex chromatographed with an apparent hydrodynamic radius of 190,000 molecular weight, the molecular weight expected for a M-CSF(rhsM-CSFR)2 complex. Other assays may be employed to measure ligand/receptor binding or receptor dimerization such as chemical crosslinking and SDS-PAGE or immunoprecipitation and SDS-PAGE. Molecules that inhibit receptor dimerization but not ligand binding provide another method to antagonize M-CSF actions.

[0064] The DNA encoding rhsM-CSFR was cloned for expression in insect cells using the following general strategy. The portion of the *c-fms* gene corresponding to amino acids one to 511 was amplified from human macrophage cDNA by polymerase chain reaction (PCR) using an upstream primer of: 5'-GCGTA<u>CCATGG</u>GCCCAGGAGTTCTGC-3' (SEQ ID NO.9) and a downstream primer of: 5'-AGTCGA<u>GGATCC</u>TCAATCCGGGGGATGCGTGTG-3' (SEQ ID NO. 10). The underlined sequences are the Ncol and BamHI restriction sites used to subclone the PCR product into the pAcC5 vector (Luckov et al, Bio/Technology 6:47-55). The pAcC5:hsM-CSFR vector was expressed in Sf9 insect cells using a baculovirus helper vector as previously described (Summers, et al., A Manual of Methods for Baculovirus Vectors and Insect Cell Culture Procedures (1987)).

[0065]Approximately two liters of serum-free 72-hour conditioned medium was collected by centrifugation and filtration from Sf9 cells infected with pAcC5:hsM-CSFR construct described above. The material was diafiltered with DEAE buffer A [10 mM Tris, pH 8.8, containing the following protease inhibitors (which were added all buffers throughout the purification): 1 mM EDTA, 2 µg/ml leupeptin and 100 µM PMSF] and concentrated 20-fold with a 20,000 molecular weight cut-off Pyrostat Ultrafiltration Membrane (Sartorius). The retentate was loaded onto a DEAE Sepharose column (Pharmacia LKB Biotechnology, Inc., Piscataway, New Jersey) having a bed volume of 100 ml that had been pre-equilibrated with DEAE Buffer A. Elution was at 5 ml/min with a 0-0.8 M NaCl gradient in 500 ml of DEAE Buffer A. Fractions enriched in rhsM-CSFR were detected by Western Analysis [Burnett, R., Anal. Biochem., 112:195 (1981)] and dot blot analysis of serially diluted fractions, using anti-c-FMS monoclonal antibodies (Oncogene Sciences, Inc.). The dot blot assay was used throughout the purification to identify fractions containing rhsM-CSF. Enriched fractions were pooled, made 0.8 M in ammonium sulfate, adjusted to pH 7.0 and loaded onto a Phenyl TSK-5-PW HPLC column (7.5 x 75 mm) (BioRad). The column was eluted at 1 ml/min with a decreasing ammonium sulfate gradient over 45 minutes, peak fractions were pooled and concentrated 10-fold with a stir cell concentrator using a YM30 membrane (Amicon). The retentate was chromatographed with FG3000XL size exclusion column (DU PONT, Wilmington, Delaware) using a mobile phase a phosphate-buffered saline (PBS) at 3 ml/min. The purified receptor was pooled, concentrated to 1 mg/ml as above and stored at 4°C. This process recovered 650 µg of rhsM-CSFR, purified 200-fold. The preparation was about 95% homogeneous as assayed by SDS-PAGE stained with Coomassie Blue.

Example 6

50

Crystallization of M-CSF/Soluble M-CSF Receptor Complex

[0066] To crystallize the M-CSF/rhsM-CSFR complexes, glass microscope slides are prepared as described in Example 2. The M-CSF composition used is incubated with a purified soluble form of the M-CSF receptor, truncated at a residue before the transmembrane region, to form an M-CSF/ receptor complex. In certain cases, the rhsM-CSFR is deglycosylated prior to the size exclusion step by incubation with N-glycanase (Genzyme, Cambridge Massachusetts)

according to the manufacturer's instructions. A small quantity of M-CSF/receptor solution is mixed in each well of the spot plate with a comparable quantity of a drop solution (such as about 20% (v/v) PEG 4000, 0.2 M MgCl₂, 0.1 M Tris-HCl (pH 8.5)). The spot plate is then placed in a clear plastic sandwich box containing a small amount of precipitant solution (such as about 23% PEG 4000, 0.2 M MgCl₂, 0.1 M Tris-HCl (pH 8.5)). The box is immediately sealed and stored at room temperature.

[0067] After a few days, crystalline M-CSF-receptor complex is isolated and redissolved in a solution containing about 50 mM Tris-HCl and is allowed to stand at room temperature. The purified M-CSF-receptor complex crystallizes from solution to form crystals for X-ray structural analysis. To facilitate solution of the crystal structure of such complexes, truncated, non-glycosylated forms of the rhsM-CSFR (described above) which retain M-CSF binding ability may be employed to generate M-CSF-receptor complex crystals.

Example 7

[0068] The biological activity of the non-glycosylated, truncated sequence used in Examples 2 and 3 was shown to be equal to that of the mature protein purified from human urine (Halenbeck, R., et al., Bio/Technology, 7:710-715 [1989]). As noted, the resulting crystals had an orthorhombic crystal lattice in the P2₁2₁2₁ space group, with cell dimensions a=33.54, b=65.26, and c=159.63 Å. Intensity data were collected using imaging plates mounted on a Weissenberg camera modified for macromolecular crystallography at the Photon Factory in Tsukuba, Japan. Native data to a nominal resolution of 2.0Å, and mercury and platinum derivative data were collected using 1.0Å radiation. Two crystal settings were used to collect native data [Rmerge (I)=7.0%, using all measurements with I>0.0].

[0069] Heavy atom derivatives of M-CSF crystals were prepared by soaking crystals in heavy atom compounds dissolved in the reservoir solution. Isomorphous and anomalous difference Patterson maps clearly revealed one site for the mercury and two sites for the platinum derivative. Anomalous and isomorphous phase information as used in initial phase refinement with the PROTEIN program package. The final figure of merit was 0.62 (8.0-3.0Å, 6960 reflections). After solvent flattening, B.C. Wang, *Methods Enzymol*, 115:90 (1985), two bundles of four alpha helices related by an approximate two-fold axis could be seen in the electron density map. Rotation and translation parameters of this noncrystallographic axis were refined by a density correlation method, J.M. Cox, *J. Mol. Biol.* 28:151 (1967). Phases were then iteratively refined by molecular avenging and solvent flattening, G. Bricogne, *Acta Cryst.*, 32:832 (1976), using an envelope calculated by putting 5Å spheres around all the atoms in the four helical bundle. Chain tracing and model building were done in the resulting map, (using the program FRODO), T.A. Jones, *Methods Enzymol.* 115:157 (1985), keeping the original MIR map as a reference.

The starting partial model for refinement contained only a polyalanine backbone for eight helices making up the two bundles. Positional refinement using the program XPLOR, A.T. Brunger, *J. Mol. Biol.* 203:803 (1985), gave an R-factor of 0.49 to 3.0Å. Phase combination with the refined MIR phases resulted in a map of sufficient quality to allow the tracing of two long loops traversing the four helical bundle and a short loop connecting two of the helices. Two strong peaks in the density, one at the top of the first helix, and the second lying directly on the molecular two fold axis, were assigned as disulfide bonded cysteines. The number of residues between these two peaks uniquely identified the position in the sequence of these cysteines and consequently the sequence of the intervening residues. This initial registration was confirmed by the presence of a number of regions of strong density corresponding to aromatic side chains in the sequence. Partial model phase combination using the added loops and those side chains that were visible allowed the remaining residues to be registered, thus determining the overall topology of the molecule. The presence of seven disulfide bonds in the dimer served as important "tether points" to confirm the correctness of the tracing.

[0071] As shown in Fig. 2, the overall topology of this form of M-CSF is that of an antiparallel four α -helical bundle, in which the helices run up-up-down-down, unlike the more commonly observed up-down-up-down connectivity of most four helical bundles. A long crossover connection links helix A to helix B and a similar connection is found between helices C and D.

[0072] A striking difference from other cytokines and other four helix bundle structures is that the truncated M-CSF α forms a disulfide-linked dimer, in which the bundles are linked end-to-end, forming an extremely flat, elongated structure (approximate dimensions 85 x 35 x 25 Å) as shown in Figs. 3 and 4A and 4B. There are three intramolecular disulfide bonds in each monomer (Cys7-Cys90, Cys48-Cys139, Cys102-Cys146) all of which are at the distal end of the molecule. One interchain disulfide bond (Cys31-Cys31) is located at the dimer interface with the noncrystallographic two-fold symmetry axis passing through it as shown in Figs. 3 and 4A and 4B. Mutation experiments indicate that all of the cysteine residues in this form of M-CSF may be necessary for full biological activity. The structure described herein suggests that their role is primarily structural rather than being related to receptor recognition.

[0073] Appendix 1 provides the three-dimensional structure of the truncated recombinant M-CSFα dimer as identified by the alpha-carbon positions of the amino acid residues in the sequence. The five carboxy terminal amino acids of each polypeptide of the dimer were not included. As will be recognized to those of skill in the art, the information in Appendix 1 is provided in the format used by the Brookhaven Protein Data Bank.

[0074] As shown, the molecule has an unusual topology which identifies important regions of M-CSF with regard to M-CSF receptor binding. Specific residues in helices A, C, and D appear to be involved in the specificity of the interaction. Altering solvent accessible residues in these regions by site directed mutagenesis to increase or decrease sidechain interactions with the receptor may be useful to generate M-CSF agonists or antagonists. For example, changing one or more histidines to non-hydrogen-donor amino acids of similar size may create an M-CSF with altered receptor binding ability.

Example 8

25

Preparation of M-CSF Heterodimers

Purification of M-CSF Monomers:

[0075] E. coli harboring the pL-M-CSF vector for NΔ3CΔ158M-CSFα described in U.S. Patent No. 4,929,700 or NΔ3CΔ221 M-CSFβ C157S, C159S (Kawasaki, et al., in Colony-Stimulating Factors, Dexter, T., Garland, J., and Testa, N., eds. [1990]) were grown in 1 liter of minimal salts medium containing glucose and the appropriate antibiotic. Expression of M-CSF was induced by shifting the temperature to 39°C for 4 hr. following the addition of casamino acids to 2%. The cells were harvested by centrifugation and lysed by sonication in 50 mM Tris (pH 8.5), 10 mM EDTA. The cell debris was recovered by centrifugation, washed with 30% sucrose containing 10 mM EDTA, and a portion of the refractile body paste was solubilized in 8 M urea under reducing conditions. After incubation at 37°C for 30 min., the solubilized M-CSF was clarified by centrifugation and filtration, and then loaded onto a Bio-Gel TSK-DEAE-5-PW column (7.5 x 75 mm) (BioRad Laboratories, Richmond, California) equilibrated in 8 M urea in 10 mM Tris (pH 8.5), 5 mM DTT, 0.5 mM EDTA. The monomeric M-CSF was eluted with a 45-min, 0-0.6 M NaCl gradient. The M-CSF peak fractions were pooled and concentrated to 10 mg/ml with a Centricon 10 microconcentrator (Amicon).

Formation and Analysis of Active M-CSF Heterodimers

[0076] The M-CSF homodimers were refolded by diluting to a protein concentration of 0.5 mg/ml in precooled 50 mM Tris (pH 8.5), 5 mM EDTA, 2 mM reduced glutathione and 1 mM oxidized glutathione, and then incubating at 4°C. The heterodimer was refolded by diluting 158 and 221F monomer pools to 1 mg/ml in the same buffer. To monitor the refolding, size exclusion high pressure liquid chromatography (SE-HPLC) analysis was performed by immediately injecting reaction samples onto a G3000SW Ultropack (Pharmacia LKB Biotechnology, Inc., Piscataway, New Jersey) column size-exclusion (7.5 x 600 mm) equilibrated in PBS (pH 6.8).

[0077] Fractionated products were analyzed on reducing SDS-PAGE and stained with Coomassie, according to the method of Laemmli, *Nature* (Canada) 227: 680-685 (1970). Biological activity was determined using the M-CSF dependent NFS-60 bioassay (see Example 10 below). Antibody neutralization experiments were carried out by pre-incubating approximately 5,000 units of M-CSF dimer with varying dilutions of the neutralizing M-CSF 5H410 Mab (made to refolded *E. coli* CΔ 150 M-CSFα dimer) prior to bioassay. (Halenbeck *et al.*, Bio/Technology 7:710 (1989)

[0078] The heterodimeric M-CSF product was designed to consist of one chain of short clone (from amino acid 4 to 158) and one chain of long clone (from amino acid 4 to 221). The long-clone chain (221F) also contained substitutions of serine for the two non-essential cysteines (at 157 and 159) to minimize the possibility of higher-order oligomer formation.

[0079] Solubilized refractile bodies of M-CSF 158 and 221F were separately chromatographed by DEAE-HPLC in 8 M urea. Only one major protein peak eluted in each case, and the peak fractions were pooled, based on an analysis of purity by non-reducing SDS-PAGE and Coomassie staining (data not shown). The resulting monomer was over 90% pure in each case. The monomers were separately concentrated to 10 mg/ml, diluted in refolding buffer, and refolded at 4°C.

[0080] To compare the rates of dimerization of short- and long-clone M-CSF, 20 µl of each refolding reaction was injected on a SE-HPLC column at 0, 2, 18 and 72 hr. The amount of dimeric M-CSF formed was determined from the peak area at the molecular weight expected for dimer. In both refolding reactions the M-CSF was mostly equilibrated to monomer at t=0 and had become about 40% dimeric by 2 hr and nearly 75% dimeric by 18 hr. The similarity of the ratio of dimer to monomer between the refolded 158 and 221F strongly suggests that the rate of dimer formation is the same for long- and short-clone M-CSF. Thus, when equal moles of 158 and 221F are present in a refolding reaction, the final relative ratios of 158 homodimer to 221F homodimer to 158/221F heterodimer are predicted to be 1:1:2 (Similar distributions have been observed *in vivo* for isozymes of lactate dehydrogenase.)

Biological Activity of Refolded Homodimers and Heterodimers

[0081] The biological activity of the refolded homodimers and heterodimers described above, was examined using the *in vitro* M-CSF-dependent NFS-60 bioassay (see Example 10 below). Figures 5A-5C show the result of these studies. These SE-HPLC and biological activity profiles analyzed after 72 hr of refolding, show that the heterodimer, Figure 5C, displays activity very similar to that of the two homodimers, (Figures 5A and 5B). Given that the separation of the heterodimer from the homodimers was nearly complete, it can be concluded that the heterodimer is fully biologically active *in vitro*.

[0082] To verify that the M-CSF protein eluting from these columns at the predicted heterodimeric position (between the two homodimers) actually did consist of equal moles of short- and long-clone monomers, analysis of a preparative purification of the 158/221F heterodimer was carried out. Phenyl-HPLC was performed as described above and was shown to completely resolve the heterodimer from the 158 and 221F homodimers, as seen in Figure 6A.

Preparative Purification of M-CSF Heterodimers:

[0083] The refolded M-CSF was adjusted to pH 7.0 with 1 N HCl, and ammonium sulfate was added to 1.2 M. The protein was loaded onto a Bio-Gel TSK-Phenyl-5-PW column (7.5 x 75 mm) (BioRad, Richmond, California) equilibrated in 1.2 M ammonium sulfate, 100 mM phosphate (pH 7.0). The M-CSF was eluted with a decreasing gradient of ammonium sulfate from 40% to 80% buffer B (10 mM phosphate, pH 7.0) in 24 min.

[0084] Reducing and non-reducing SDS-PAGE (Figures 6B and C) showed that internal controls (the 158 and 221F dimers) were purified to approximately 95% homogeneity by this column, and each consisted of the single expected monomeric band. The gel analysis also showed that the heterodimer was purified to approximately 95% homogeneity and that it consisted of equivalent amounts of 158 and 221F monomers. Recovery of purified 158/221F heterodimer from refractile body paste to final product was greater than 15 %.

²⁵ [0085] The bioactivity of the dimeric M-CSF species was determined and, when compared to the A₂₈₀ profile in Figure 6A, confirms the finding that the heterodimer is fully active. The specific activity of the 158/221F heterodimer, calculated using the peak fraction, was 8.0 x 10⁷ units/mg, compared to 9.0 x 10⁷ and 6.8 x 10⁷ units/mg for 158 and 221F homodimers, respectively.

[0086] The biological activity of all three dimer species was neutralized to the same extent in serial dilution neutralization experiments using the 5H410 M-CSF Mab in the NFS-60 bioassay. This antibody also neutralizes "naturally refolded" Chinese hamster ovary cell (CHO)-expressed M-CSF in a similar fashion. This result further suggests that the refolded conformation of the new M-CSF heterodimer is essentially native-like, at least with regard to the region within the first 150 amino acids that is responsible for *in vitro* activity.

35 Example 9

15

Selection of Amino Acid Substitutions in M-CSF Based on Crystallographic Data

[0087] The X-ray crystallographic data described above provided sufficient structural information regarding M-CSF to be able to identify a limited subset of the amino acids in the protein that are likely to be crucial for M-CSF receptor binding and biological activity and thus which represented likely candidates for mutagenesis with the ultimate goal of providing M-CSF muteins having altered biological activity (*i.e.*, agonists or antagonists). Based on this information, several criteria were used to generate a list of possible target amino acids for substitution.

[0088] The first criterion was solvent exposure or solvent accessibility, which refers to amino acids residues at the surface of the protein. Residues having a solvent accessible surface area of greater than about 0.25 and preferably greater than about 0.4 are preferred based on normalization of the surface area of the amino acid accessible when in the trypeptide gly-x-gly (Kabsch, W. et al., Biopolymers 22:2577 (1983)). Residues were chosen which do not interact with other parts of the protein such as the dimer interface in order to maintain the relative orientation of monomers and to avoid disturbing the process of protein folding. Still another criterion used in certain instances in selecting candidate amino acid substitutions is the relationship of the residues to corresponding residues in mouse M-CSF. Another important selection criterion was that the substitutions be non-conservative so as to attempt to disrupt possible hydrogen bonding or hydrophobic interactions with M-CSF receptor residues.

[0089] Table 1 lists exemplary amino acid residues and exemplary substitutions. Using the criteria for selecting candidates for substitutions set forth above, those of ordinary skill in the art may readily ascertain other possible candidates for substitution.

Table I

Table 1						
	Candidate Substitutions					
Wild Type An	nino Acid and ation	Substitutions				
His (H)	15	Ala(A) or Leu(L)				
Gin (Q)	17	Ala(A) or Glu(E)				
Gln (Q)	79	Ala(A) or Asp(D)				
Arg (R)	, 86	Glu(E) or Asp(D)				
Glu (E)	115	Ala(A)				
Glu (E)	41	Lys(K) or Arg(R)				
Lys (K)	93	Ala(A) or Glu(E)				
Asp (D)	99 -	Lys(K) or Arg(R)				
Leu (L)	55	Gin(Q) or Asp(N)				
Ser (S)	18	Ala(A) or Lys(K)				
Gln (Q)	20	Ala(A) or Asp(D)				
Arg (R)	21	Ala(A), Glu(E), or Asp(D)				
lle (i)	75	Lys(K) or Glu(E)				
Val (V)	78	Lys(K) or Arg(R)				
Leu (L)	85	Glu(E) or Asn(N)				
Asp (D)	69	Lys(K) or Arg(R)				
Asn (N)	70	Ala(A) or Glu(E)				
His (H)	9	Ala(A) or Asp(D)				
Asn (N)	63	Lys(K) or Arg(R)				
Thr (T)	34	Gln(Q) or Lys(K)				

[0090] It is not expected that every candidate substitution listed will result in the production of M-CSF agonists or antagonists (see Example 12 below). Rather they represent a non-exclusive list of candidates likely to result in the production of agonists or antagonists based on the selection criteria set forth above. It should also be noted that even if a variant does not act as an agonist or antagonist when compared with native M-CSF, the variant is still useful for conventional uses of the ligand (if it retains the same activity as the ligand) or as for example, a diagnostic reagent.

Example 10

10

15

20

25

35

45

Preparation of H9A, H15A M-CSF Muteins

[0091] A variety of M-CSF muteins with altered solvent-accessible residues from regions of the M-CSF mature N terminus and helices A, C, and D were constructed using techniques known in the an. For example, two histidines in the N-terminal/A helix region were changed to alanine through site-directed mutagenesis of a truncated form of M-CSFα (encoded by pLCSF158A). Involvement of one of three M-CSF histidine residues in M-CSF receptor interaction was implicated by our observation that diethylpyrocarbonate (DEPC) modification of histidines in M-CSF at a 1:100 DEPC:histidine ratio (as described in *Meth. in Enzymol.* 47:431 (1977)) significantly reduced bioactivity.

[0092] Plasmid DNA pLCSF158A was prepared from the *E. coli* strain HW22 carrying the plasmid pLCSF158A (U.S. Patent No. 4,929,700, Example 6, "*E. coli* strain HW22 transformed with pJN653 containing the asp₅₉SCSF/NΔ3CΔ158 gene"). The strain was grown in 350 ml R2 media (2X Luria Broth containing 1% sodium chloride and no glucose, *J. Bact.*, 74:461 (1957)) containing 50 micrograms/ml ampicillin at 30° C with shaking overnight. Plasmid DNA was prepared from the cells using a Qiagen-tip 100 column according to the manufacturer's directions.

[0093] Twenty micrograms of pLCSF158A DNA were digested with 66 units of *Hind*III and 66 units of *Stu*I at 37°C for 3hr. 20 min. in 200 microliters 1X New England Biolabs NEBuffer 2 (New England Biolabs, Beverly, Massachusetts). The DNA was extracted with phenol and chloroform, then ethanol precipitated. The DNA was treated with one unit of Calf Intestinal Alkaline Phosphatase in 100 microliters of 1X Dephosphorylation Buffer, supplied by Boehringer Mannheim (Indianapolis, Indiana), at 37°C for 30 min. An additional unit of Calf Intestinal Alkaline Phosphatase was added to the reaction and incubation was continued at 50°C for 1 hr.

[0094] The resulting DNA was then run on a 1 % FMC Bioproducts (Rockland, Maine) Sea KEM® GTC® agarose gel. The 5.7 kb pLCSF158A fragment was cut from the gel and purified on Qiagen (Chatsworth, California) Qiaex beads according to the manufacturer's directions.

[0095] Polymerase chain reaction (PCR) was then performed and a PCR product was produced that contained a mutagenized M-CSF sequence in which histidines 9 and 15 (counting from the mature N-terminus) were altered to alanine (generating an H9A, H15A PCR fragment). The 5' portion of the M-CSF gene was amplified from the plasmid pLCSF158A in a PCR reaction using the primers LF73 and LF74. Details of PCR are provided by Mullis, K. et al., U.S. Patent No. 4,683,202; Ehrlich, H., U.S. Patent No. 4,582,788; Saiki et al., U.S. Patent No. 4,683,195; Mullis, K.B., Cold Spring Harbor Symp. Quant. Biol. 51:263-273 (1986); Saiki et al., Bio/Technology 3:1008-10012 (1985); and Mullis, K.B. et al., Meth. Enzymol 155:335-350(1987), all of which are incorporated herein by reference. The sequences of these primers are:

LF73:

5'-AGGTGTCTCATAGAAAGTTCGGACGCAGGCCTTGTCATGCTCTTCATAATCCTTGG-3' (SEQ ID NO. 1)

LF74:

5'-CAGGAGAAAGCTTATGTCTGAATATTGTAGCGCCATGATTGGGAGTGGAGCCCTGCAG-3' (SEQ ID NO. 2)

30

25

20

[0096] The expected PCR product was designed to include 337 bp of pLCSF158A sequence, with the *Hind*III and *Stul* sites located at each end of the product for cloning, and the histidine codons for His 9 and His 15, CAC, mutated to an alanine codons, GCC.

[0097] This product was amplified in four separate PCR reactions each containing 100 ng of pLCSF158A DNA, 50 pmoles LF73; 50 pmoles LF74, 37.5 μM dNTPS, 5% glycerol, 1X Perkin Elmer Cetus PCR Buffer, and 2.5 units of Perkin Elmer Cetus AmpliTaq[®] DNA Polymerase in a 100 microliter volume. The amplification was carried out in a Perkin Elmer Cetus DNA thermocycler. Before adding the AmpliTaq[®], the reactions were brought to 95°C. The amplification was carried out for 25 cycles ramping to a denaturation temperature of 95°C in 1 sec., denaturing at 95°C for 1 min;

ramping to an annealing temperature of 68°C in 1 sec., annealing at 68°C 1 min.; ramping to an extension temperature of 72°C in 30 sec., extending at 72°C 1 min., 30 sec. Final extension was carried out at 72°C for 10 min.

[0098] Five microliters of each reaction were run on a 3% agarose gel (1.5% FMC Bioproducts SeaKem[®] GTG[®] agarose, 1.5% FMC Bioproducts NuSeive[®] GTG[®] agarose in Tris-Borate buffer) (FMC Bioproducts, Rockland, Maine). Gels were then stained with ethidium bromide. For each reaction, a major band was visible at approximately 337 bp.

- The four reactions were pooled, extracted with phenol and chloroform, precipitated with ethanol, resuspended and digested with 250 units of *Stul* in a final volume of 500 microliters 1X NEBuffer 2 at 37°C for 2hr., 500 units of *Hind*III were added to the reaction, the volume increased to 1 ml in 1X NEBuffer 2 and digestion was continued at 37°C for an additional 2.5hr. The DNA was electrophoresed on a 3% agarose gel. The 300 bp digested product was cut from the gel and purified on Qiagen Qiaex beads according to the manufacturer's directions.
- 50 [0100] Approximately 68 ng of the *Hind*III/Stul digested PCR product was then ligated to approximately 28 ng of the 5.7 kb *Hind*III/Stul digested pLCSF158A vector DNA at an insert-to-vector ratio of approximately 5:1. Ligation was carried out with 1 unit of Boehringer Mannheim T4 DNA ligase in 1X ligation buffer, supplied by the manufacturer, in a 20-microliter volume at 16°C overnight. As a control 28 ng of the 5.7 kb *Hind*III/Stul digested pLCSF158A vector DNA was ligated to itself under the same conditions with no insert present.
- 55 [0101] Half of each ligation mixture was used to transform competent E. coli DG116 (ATCC# 53606) cells using a protocol similar to the calcium chloride procedure described in Molecular Cloning a Laboratory Manual Maniatis et al., Cold Spring Harbor Laboratory (1982). Transformed cells were allowed to express at 30°C with no selection for 90 min., plated on R2-4 (10g tryptone, 5g yeast extract, 5g NaCl, 2 drops antifoam A, 4 ml 50% glucose and 15g agar in 1 liter).

plates containing 50 micrograms/ml ampicillin. The plates were incubated at room temperature 72 hr. One fourth of each transformation was plated. For the ligation containing the insert, 66 ampicillin resistant colonies appeared on the plates. For the ligation with no insert, no colonies appeared.

[0102] One of these colonies, designated strain TAF172-2, was picked and cultured in 350 ml R2 broth with 50 micrograms/ml ampicillin at 30°C with shaking overnight. A frozen stock in 40% glycerol was made from this culture and stored at -70°C. DNA was isolated from the culture using Qiagen-tip 100 columns as described above.

[0103] The purified DNA, pTAF172-2, was sequenced using the di-deoxy method and shown to contain the sequence of pLCSF158A coding for M-CSF N∆3C∆158 with His 9 and His 15 mutated to encode Ala.

[0104] The M-CSF mutein NΔ3CΔ158 H9A, H15A encoded by pTAF172-2 was expressed, purified, refolded to form dimeric protein and assayed essentially as described in U.S. Patent No. 4,929,700 Example 5, using 8M urea as a denaturant and in the DEAE purification step.

[0105] The N-terminal sequence of the purified mutein was determined through 20 cycles, using a standard automated Edman degradation method, and shown to be identical to that of the parental N Δ 3C Δ 158 M-CSF α reference protein except that His 9 and His 15 had been altered to Ala. Protein concentration was determined using A280 and an extinction co-efficient of 0.68.

[0106] The purified mutein dimers were subjected to bioassay using NFS-60 cells which is an M-CSF dependent cell line which forms colonies in the presence of active M-CSF. Standards and purified mutein samples were serially diluted 1:2 in RPMI media (10% fetal bovine serum) in a 50 microliter volume in 96-well microtiter plates. 50 microliters NFS-60 cells (ATCC NO. CRL 1838), maintained in 4000U/ml M-CSF, washed 2X, and diluted to a concentration of 1 x 10⁵ cell/ml, were added to each sample well. Plates were incubated at 37°C, 5% CO₂ for 44h. 5 mg/ml MTT dye was added to each sample and incubation continued for 4h. 20% SDS was added to each sample. The plates were wrapped in foil and left overnight. The absorbance of each sample was read at 570 nm and compared to a standard known M-CSF concentration. The H9A, H15A mutein showed a specific activity 7.6x10³ U/mg compared to 6.9x10⁷ U/mg for the parental M-CSF NΔ3CΔ158 reference in the same assay. This represents a nearly 10,000 fold reduction biological activity for the mutein. The same M-CSF mutein preparation was shown to have greatly decreased M-CSF receptor-binding ability using the NFS-60 receptor competition assay described in this example. Because the H9A, H15A M-CSF mutein was otherwise not significantly different from the parental M-CSFα, including crystallizability and space group (see Example 12) we believe the decrease in biological activity is not due to gross deformation of structure but reflects an alteration in important M-CSF receptor contacts.

[0107] Using essentially the same methodology, two M-CSF muteins contacting singly substituted histidines at residues 9 and 15 were generated (e.g. H9A and H15A). The H9A construct utilized LF80:

5'-CAGGAGAAAGCTTATGTCTGAATATTGTAGCGCCATGATTGGGAGTGGACACCTGCAG-3' (SEQ ID NO 3);

and LF73 (described in this example) as the PCR primers. The H15A construct utilized LF81:

35

40

5'-CAGGAGAAAGCTTATGTCTGAATATTGTAGCCACATGATTGGGAGTGGAGCCCTGCAG-3'
(SEQ ID NO. 4);

and LF 73 (described in this example) as PCR primers. Biological assay of the purified muteins immediately following the refolding step described above showed approximate biological specific activities as follows: 4 x 10⁶ U/mg for H9A and less than 3 x 10³ U/mg for H15A, in an assay which the parental M-CSF construct displayed 8 x 10⁷ U/mg. This information, combined with that described above, suggests that H15A as well as possible H9A represents contacts that are important for M-CSF receptor binding. Nearby solvent accessible residues such as Y6 and S13 (see also Table 1) may also represent M-CSF receptor contact residues. Non-proteinaceous mimics of the side chains of H9, H15, and nearby solvent accessible side chains may represent M-CSF agonists or antagonists. Such residues should be left unchanged in M-CSF mutein constructs designed to retain full M-CSF receptor binding but to have M-CSF antagonist properties because they lack significant M-CSF bioactivity. Homodimers of muteins that retain full receptor-binding ability and display significantly reduced bioactivity should represent M-CSF antagonists. M-CSF muteins that are greatly decreased in both M-CSF bioactivity and receptor binding ability (such as H15A) may generally be useful in M-CSF immunoassay applications and might represent useful therapeutic agents for patients having auto-antibodies to M-CSF immunoassay applications and might represent useful therapeutic agents for patients having auto-antibodies to M-CSF.

Example 11

Preparation of Q20A, V78K M-CSF Muteins

[0108] Using essentially the same methodology described in Example 10, a double mutant of M-CSF (Q20A, V78K) was constructed to test the importance of solvent accessible residues in the central portion of helices A and C. The following PCR primers were used.

LF63:

10

15

20

5'-AGGAGAAAGCTTATGTCTGAATATTGTAGCCACATGATTGGGAGTGGACACCTGCAGT CTCTGGCTCGGCTG-3' (SEQ ID NO. 5)

LF64:

5'-GGACGCAGGCCTTGTCATGCTCTTCATAATCCTTGGTGAAGCAGCTCTTCAGCCTCAA AGAGAGTTCCTGCAGCTGTTTAATGGC-3' (SEQ ID NO. 6)

[0109] The resulting mutein was expressed, refolded, purified, and assayed as described in Example 8. The specific biological activity was 1.4 X 10^7 U/mg, approximately 8-10 fold lower than that of the parental M-CSF α reference standard. The receptor binding activity of this mutein was also decreased.

[0110] This result again supports the prediction of over crystallographic study of truncated M-CSFα which concluded that important M-CSF receptor contact residues exist among the solvent accessible residues in helices A and/or C and/or D. Certain of these mutations will, as we have shown, have lower biological activity and lower M-CSF receptor-binding ability. Some may have lower biological activity without a decrease in receptor-binding ability. Some may have increased biological activity and receptor binding ability, and some may have no affect on either.

[0111] Two examples of the latter are Q17A, R21A (produced using PCR primers LF72:

5'-TTGTAGCCACATGATTGGGAGTGGACACCTGGCGTCTCTGCAGGCGCTGATTGAC-3'

and LF73 (described in Example 9) and E115A, N119A (produced using LF75:

5'-CATGACAAGGCCTGCGTCCGAACTTTCTATGAGACACCTCTCCAGTTGCTGGCGAAGG TCAAGGCTGTCTTTAATG-3' (SEQ ID NO. 7);

and

40

45

LF79: 5'-GGATCAGGATCCCTCGGACTGCCTCTC-3' (SEQ ID NO. 8)).

Both of these constructs changed side chain properties of solvent-accessible amino acids in the areas of interest but did not affect biological specific activity, compared to the parental reference molecule. These results indicate that residues Q17, R21, E115, and N119 do not need to be altered in muteins designed to have M-CSF agonist or antagonist activity. In fact, to minimize the likelihood of antibody formation to potentially administered M-CSF-based proteinaceous drugs, it is desirable to retain the solvent-accessible parental M-CSF residues (to resemble the native molecule) whenever possible.

[0112] The retained activity of the muteins including changes at Q17, R21, E115, and N119 does not rule out large effects on activity contributed by nearby residues (such as H15). In fact, the regions we have altered are predicted by the crystal structure to be important for receptor binding and/or signaling. Antagonistic M-CSF muteins may require use of multiple residue changes per mutein or use of heterodimeric molecules containing one or more mutations in each polypeptide chain, since M-CSF residues important in receptor signaling are believed to be composed of discontinuous regions of M-CSF.

Example 12

Formation of M-CSF Heterodimers Having Decreased Receptor-Binding Ability and/or Decreased Biological Specific Activity

[0113] M-CSF can be folded *in vitro* to generate fully active heterodimers, as shown in Example 8. By making heterodimers of M-CSF which incorporate M-CSF muteins with altered M-CSF signaling ability, it should be possible to generate antagonists of M-CSF useful for treatment of patients with M-CSF mediated diseases. To generate a heterodimer containing one subunit of M-CSFαNΔ3/CΔ158, H9A/H15A and one subunit of M-CSFβNΔ3/CΔ221, C157S/C159S, each mutein was expressed in *E. coli* and purified separately by DEAE Sepharose under denaturing and reducing conditions as described in Example 8. The two muteins subunits were mixed together prior to refolding to generate a solution containing a final mutein molar ratio of 1:1, then this solution was diluted to 0.2 mg/ml with refolding buffer as described in Example 8. Following refolding, the heterodimeric molecule was separated from the homodimers by two consecutive passes over a Phenyl TSK-5-PW HPLC column as described in Example 8. No contaminants were detected when the purified heterodimer preparation was examined, by non-reduced SDS-PAGE or size exclusion HPLC using a BioSil SEC250 column (BioRad).

The purified heterodimer was submitted to the NFS60 cell based bioassay described in Example 8. The cal-[0114] culated specific activity was 2.9 x 10⁶ U/mg which correlated to a 35-fold reduction as compared to the activity of the parental M-CSF heterodimer described in Example 10. The relative binding affinity to cell surface M-CSF receptor was measured by radioligand displacement in which the displacement of 125I-M-CSF from an M-CSF receptor by an M-CSF mutein was measured using methods well known in the art. In brief, the following were added in a final volume of 100µl in each well of a 96-well cell culture plate; approximately 80,000 cpm of purified recombinant human M-CSF labeled with 1251 (using lodobeads as described by the manufacturer, Pierce, Rockford, IL), 300,000 NFS-60 cells that had been washed and then grown for 18 hours in growth medium minus the normal maintenance level of M-CSF, plus unlabeled M-CSF that had been serially diluted in the same medium. The plates were incubated at 4°C for 20 hours and the cells were collected on glass-fiber filters. Maximum binding was measured in the absence of unlabeled M-CSF and non-specific binding was measured in the presence of 1000-fold greater concentration of unlabeled M-CSF (compared to labeled M-CSF). The concentration of M-CSF required to inhibit 50% of the labeled M-CSF binding to the cells (IC50) was used to determine differences in affinity. Results are expressed as percent displacement of radioactive M-CSF versus mutein concentration (Figure 7). The IC_{50} of the heterodimer (Fig. 7 closed squares) was reduced 30-fold to about 500 pM as compared to an IC₅₀ of about 17 pM for M-CSFαNΔ 3/CΔ158 (158) (Fig. 7 closed circles). The similarity between the reduction in specific activity and receptor affinity of the heterodimer indicates that the reduction in bioactivity was due to decreased receptor-binding ability. Similarly, the binding affinities of the Q20A, V78KF (Fig. 7 open circles) and H9A/H15A (open squares) muteins were also measured in this radioligand displacement assay. The Q20A,V78K mutein had an IC50 of about 100 pM and the H9A,H15A mutein has an IC50 of about 1 µM, correlating to decreased binding affinities of 5-fold and 50,000-fold, respectively. For each mutein, the reduction in receptor affinity was similar to the reduction in specific activity, again indicating that the reduction in bioactivity was due to reduced receptor-binding ability.

40 Example 13

Crystallization and Characterization Of M-CSF H9A, H15A Muteins

[0115] The H98, H15A mutein described in Example 10 was crystallized using the hanging drop method described in Examples 1 and 2 using the following buffer conditions: 30% polyethylene glycol 4000; 100 mM Li $_2$ SO $_4$; and 100 mM Tris pH 8.5. The crystals produced under these conditions were mombohedral prisms having dimensions of 0.7 mm x 0.2 mm x 0.2 mm. X-ray crystallographic analysis using precession photographs showed crystals in the P2 $_1$ 2 $_1$ 2 $_1$ 1 space group with cell dimensions of a = 33.99, b = 65.37, c = 159.90, d = 90, e = 90, and f = 90 angstroms and diffract to a nominal resolution of 3 angstroms. These physical properties are essentially the same as those observed for the parental NA3C Δ 158 M-CSF α molecule and suggests that the biological effects of the H9A, H15A alterations are not the consequence of gross global alterations in M-CSF structure, but rather are the result of altered side chains that are important in interacting with the M-CSF receptor. Alteration of those histidine side chains may have affected receptor binding by changing atoms that interact with, stabilize or facilitate receptor binding or changes in receptor conformation. Changes such as H15A may also have affected these functions by altering the position of the nearby side chain in M-CSF, most likely in the A and/or C helix regions.

[0116] The foregoing examples are presented by way of example and are not intended to limit the scope of the invention as set forth in the appended claims.

APPENDIX 1

10	MOTA	10 CA SER	4	63.753 80.590 222.385 1.00 58.89
10	MOTA	13 CA GLU	5	64.883 77.664 219.972 1.00 59.56
	MOTA	23 CA TYR	6	62.285 76.840 217.324 1.00 54.54
	MOTA	37 CA CYS	7 .	63.509 80.109 215.834 1.00 54.96
15	MOTA	44 CA SER	8	66.853 79.529 214.160 1.00 54.41
	MOTA	52 CA HIS	9	65.466 77.798 211.053 1.00 55.11
	MOTA	65 CA MET	10	61.857 78.767 211.073 1.00 50.66
	MOTA	74 CA ILE	11	62.173 80.905 207.970 1.00 48.31
20	ATOM	83 CA GLY	12	63.487 78.618 205.354 1.00 52.83
	ATOM	88 CA SER	13	64.609 79.967 201.952 1.00 53.13
	MOTA	96 CA GLY	14	61.665 78.009 200.514 1.00 46.05
•	MOTA	101 CA HIS	15	59.455 80.924 201.933 1.00 40.25
25	MOTA	114 CA LEU	16	62.210 83.200 200.992 1.00 41.58
	MOTA	123 CA GLN	17	62.153 82.266 197.328 1.00 46.17
	MOTA	135 CA SER	18	58.378 81.868 197.227 1.00 46.55
30	ATOM	143 CA LEU	19	58.280 85.544 198.199 1.00 44.75
	MOTA	152 CA GLN	20	60.814 86.044 195.458 1.00 39.32
	MOTA	164 CA ARG	21	58.457 84.718 192.736 1.00 36.91
	MOTA	181 CA LEU	22	55.610 86.694 194.252 1.00 38.27
35	MOTA	190 CA ILE	23	57.676 89.822 193.465 1.00 34.84
	MOTA	199 CA ASP	24	58.677 88.274 190.086 1.00 31.50
	MOTA	208 CA SER	25	55.086 87.802 188.978 1.00 31.73
	MOTA	216 CA GLN	26	54.284 91.540 189.154 1.00 33.24
40	MOTA	228 CA MET	27	53.961 93.634 186.082 1.00 34.82
	MOTA	237 CA GLU	28	56.227 96.566 186.089 1.00 36.54
	MOTA	247 CA THR	29	53.758 99.303 186.920 1.00 44.80
45	ATOM	256 CA SER	30	54.590 102.747 188.323 1.00 49.40
40	MOTA	264 CA CYS	31	51.427 103.423 190.046 1.00 43.03
	MOTA	271. CA GLN	32	52.166 104.336 193.613 1.00 33.74
	MOTA	283 CA ILE	33	50.430 102.961 196.634 1.00 32.48
50	MOTA	292 CA THR	34	50.240 104.213 200.133 1.00 38.94
	ATOM	301 CA PHE	35	51.291 102.173 203.199 1.00 39.16
	MOTA	313 CA GLU	36	52.707 102.849 206.761 1.00 32.09
-	MOTA	323 CA PHE	37	56.073 101.358 207.299 1.00 25.78

20

```
59.082 101.751 209.490
                                                               1.00 34.46
                              38
      MOTA
              335 CA VAL
                                                               1.00 44.03
                                      61.044 104.798 208.714
              343
                  CA ASP
                              39
      MOTA
5
                                                               1.00 54.20
                  CA GLN
                              40
                                      64.648 103.690 208.314
      MOTA
              352
                                      65.924 107.142 209.332
                                                               1.00 52.01
              364
                  CA GLU
                              41
      MOTA
                                      63.934 107.629 212.631
                                                               1.00 44.26
                  CA GLN
                              42
      ATOM
              374
                                      64.770 104.161 213.955
                                                               1.00 45.54
                  CA LEU
                              43
              386
      MOTA
10
                                      68.126 102.952 212.789
                                                               1.00 51.53
                   CA ALA
                              44
      MOTA
              395
                                      69.175 100.197 215.232
                                                               1.00 47.92
              401
                   CA ASP
                              45
      MOTA
                                                               1.00 45.27
                                              97.054 213.098
                                      68.861
      MOTA
              410
                   CA PRO
                              46
                                              94.587 215.613
                                                               1.00 39.92
                   CA VAL
                              47
                                      67.352
15
      MOTA
              417
                                      64.862 .97.220 216.692
                                                               1.00 36.44
                   CA CYS
                              48
      MOTA
              425
                                                               1.00 37.39
                                             98.158 213.105
              432 CA TYR
                              49
                                      64.089
      MOTA
                                              94.474 212.189
                                                               1.00 36.15
                   CA LEU
                              50
                                      63.436
      MOTA
              446
20
                                              94.029 215.472
                                                               1.00 38.66
                                      61.603
              455
                   CA LYS
                              51
      MOTA
                                                               1.00 41.72
                              52
                                      59.322
                                              96.862 214.367
      MOTA
              468
                   CA LYS
                                              96.177 210.534
                                                               1.00 35.81
                   CA ALA
                              53
                                      59.229
      MOTA
              481
                                                               1.00 40.82
                                              92.701 211.432
                                      58.096
                   CA PHE
                               54
      MOTA
              487
25
                                              93.733 213.614
                                                               1.00 42.59
                                      55.236
      MOTA
              499
                   CA LEU
                              55
                                              95.869 210.697
                                                               1.00 44.24
      MOTA
              508
                   CA LEU
                               56
                                      54.205
                                              93.257 207.992
                                                               1.00 35.10
                                      54.596
                   CA VAL
                               57
      MOTA
              517
30
                                              91.431 209.628
                                                               1.00 39.91
                                      51.646
              525
                   CA GLN
                               58
       MOTA
                                               94.254 208.663
                                                               1.00 39.60
                                      49.277
       MOTA
              537
                   CA ASP
                               59
                                                               1.00 37.05
                                              94.729 205.185
                                      50.755
       MDTA
               546
                    CA ILE
                               60
                                                               1.00 33.86
                                      50.303
                                               90.932 204.676
                    CA MET
                               61
       MDTA
               555
35
                                               90.931 206.028
                                                                1.00 48.10
                    CA GLU
                                       46.766
                               62
       MDTA
               564
                                               94.198 204.286
                                                                1.00 49.51
                    CA ASP
                                       45.602
                               63
       ATOM
               574
                                               95.244 201.187
                                                                1.00 42.65
                       THR
                               64
                                       47.659
       MOTA
               583
                    CA
                                               92.155 199.603
                                                                1.00 37.90
40
                    CA MET
                               65
                                       49.106
       MOTA
               592
                                               91.249 198.376
                                                                1.00 39.50
                    CA ARG
                                       45:646
       MOTA
               601
                               66
                                                                1.00.35.13
                                               88.245 195.824
                    CA PHE
                               67
                                       45.516
       MDTA
               618
                                               86.013 194.326
                                                                1.00 44.68
       MOTA
               630
                    CA ARG
                               68
                                       42.655
45
                                               83.454 196.865
                                                                1.00 47.42
                    CA ASP
                                       41.727
       ATOM
               647
                               69
                                               80.057 195.984
                                                                1.00 44.47
                    CA ASN
                               70
                                       43.117
       MOTA
               656
                                                                1.00 34.80
                                       45.873 81.088 193.518
       MOTA
               667
                    CA THR
                               71
```

55

```
MOTA
             676 CA PRO
                            72
                                    49.607
                                           80.251 194.114 1.00 29.78
      MOTA
             683 CA ASN
                            73
                                    50.307
                                           83.745 195.418
                                                           1.00 28.47
      MOTA
             694 CA ALA
5 .
                            74
                                    47.469
                                           83.922 197.953
                                                           1.00 23.24
      MOTA
             700 CA ILE
                            75
                                           80.431 199.144 1.00 22.87
                                    48.690
      MOTA
             709
                  CA ALA
                            76
                                    52.212
                                           81.892 199.669
                                                           1.00 29.47
      MOTA
             715 CA ILE
                            77
                                    50.816
                                           84.796 201.668
                                                           1.00 37.21
10
      MOTA
             724
                 CA VAL
                            78
                                    48.902
                                           82.158 203.857
                                                           1.00 33.69
      MDTA
             732 CA GLN
                            79
                                    52.207
                                           80.382 204.466 1.00 33.36
      MOTA
             744
                 CA LEU
                            80
                                   54.040 83.518 205.449 1.00 32.55
      MOTA
             753
                 CA GLN
                                   51.152 84.517 207.743 1.00 35.07
                            81
15
      MOTA
             765 CA GLU
                            82
                                    51.541
                                           81.115 209.443 1.00 37.26
      MOTA
                 CA LEÚ
             775
                            83
                                    55.367 81.348 209.619
                                                           1.00 35.10
      MOTA
             784
                 CA SER
                            84
                                           84.840 211.066 1.00 39.64
                                    54.960
20
      MOTA
             792
                 CA LEU
                            85
                                    52.608 83.541 213.704
                                                           1.00 42.55
      MOTA
             801
                  CA ARG
                            86
                                    54.984 80.903 215.121 1.00 42.96
      MOTA
             818
                 CA LEU
                            87
                                    57.916 83.317 214.670
                                                           1.00 45.11
      MOTA
             827
                 CA LYS
                            88
                                    56.504 85.876 217.138
                                                           1.00 50.71
25
      MOTA
             840 CA SER
                            89
                                    57.621 83.610 219.825
                                                           1.00 54.58
      MOTA
             848
                  CA CYS
                            90
                                    61.211 84.524 219.014
                                                           1.00 46.79
      MDTA
             855
                 CA PHE
                            91
                                    60.235
                                           88.070 219.743 1.00 46.37
     MDTA
             867 CA THIR
                            92
                                   59.619 89.812 223.026 1.00 52.01
30
      ATOM:
             876
                 CA ALA
                            93
                                    56.916 92.445 223.523
                                                           1.00 56.09
      MOTA
             882
                 CA ASP
                            94
                                    57.914 96.123 223.998 1.00 59.35
      MOTA
             891 CA TYR
                            95
                                    55.685 99.100 225.194
                                                           1.00 66.22
35
      MOTA
             905
                  CA GLU
                            96
                                    52.401 100.066 223.612 1.00 65.51
      MOTA
             915 CA GLU
                            97
                                    53.343 103.890 223.182
                                                           1.00 62.60
      MOTA
             925
                 CA HIS
                            98
                                  56.046 102.712 220.836 1.00 58.47
      MOTA
             938
                 CA ASP
                            99
                                    53.422 100.733 218.887 1.00 59.97
      MDTA
             947
                 CA LYS
                           100
                                    52.950 103.470 216.162 1.00 57.00
      MOTA
             960
                 CA ALA
                           101
                                    56.259 105.270 216.975 1.00 48.82
      MDTA
             966 CA CYS
                           102
                                    58.861 105.450 214.154 1.00 42.41
      MOTA
             973
                 CA VAL
                           103
45
                                    56.175 104.773 211.621
                                                           1.00 33.87
      MOTA
             981
                 CA ARG
                           104
                                    55.993 106.851 208.456 1.00 42.76
      MOTA
             998 CA THR
                           105
                                   53.577 106.853 205.501 1.00 38.00
```

50

```
MOTA
           1007 CA PHE
                           106
                                   54.931 105.944 202.136
                                                           1.00 33.48
     MOTA
           1019
                 CA TYR
                           107
                                   53.853 106.546 198.455
                                                           1.00 34.28
5
     MOTA
           1033
                CA GLU
                           108
                                   55.765 103.906 196.543
                                                           1.00 33.88
    MOTA
           1043
                 CA THR
                           109
                                   55.266 101.309 193.854
                                                            1.00 36.94
    MOTA
           1052 CA PRO
                           110
                                   53.745 97.857 194.657
                                                            1.00 29.84
     MOTA
           1059 CA LEU
                           111
                                   57.132 96.374 193.761
                                                            1.00 35.66
10
     MOTA
           1068 CA GLN
                           112
                                   58.991 98.756 196.160
                                                           1.00 33.48
     MDTA
           1080 CA LEU
                                          97.691 198.912
                           113
                                   56.612
                                                            1.00 24.26
            1089 CA LEU
                                   57.253 93.932 198.091
     MOTA
                           114.
                                                            1.00 30.75
15
     MOTA
            1098 CA GLU
                           115
                                   61.063
                                          94.385 197.984
                                                            1.00 39.68
     MOTA
            1108 CA LYS
                           116
                                   60.805 95.824 201.526
                                                            1.00 39.24
     MOTA
           1121 CA VAL
                           117
                                   58.669 92.840 202.743
                                                            1.00 32.46
     MOTA
            1129
                 CA LYS
                           118
                                   61.312 90.620 201.111
                                                            1.00 37.34
20
     MOTA
            1142
                 CA ASN
                           119
                                   64.231 92.316 202.959
                                                            1.00 39.46
     MOTA
                                   62.470 92.033 206.370
            1153 CA VAL
                           120
                                                            1.00 32.07
     MOTA
                                   61.935 88.268 206.034
            1161 CA PHE
                           121
                                                            1.00 26.35
            1173 CA ASN
                                    65.404 88.049 204.416
     MOTA
                           122
                                                            1.00 36.04
25
                                    67.372 89.816 207.138
     MOTA
            1184 CA GLU
                           123
                                                            1.00 40.82
     MOTA
            1194 CA THR
                           124
                                           88.435 210.122
                                    65.331
                                                            1.00 41.72
                                    66.456 85.000 208.648
                           125
     MOTA
            1203 CA LYS
                                                            1.00 43.07
30
                                    70.010 86.339 208.179
     MOTA
            1216 CA ASN
                           126
                                                            1.00 48.89
     MDTA
            1227 CA LEU
                           127
                                    70.000 87.300 211.903
                                                            1.00 47.97
     MOTA
                           128
                                            84.194 213.435
                                                            1.00 47.21
            1236 CA LEU
                                    68.315
     MDTA
                           129
                                            82.134 211.864
                                                            1.00 50.11
            1245 CA ASP
                                    71.086
35
     MOTA
                                            84.146 214.117
                                                            1.00 52.04
            1254 CA LYS
                           130
                                    73.448
     MOTA
                                            84.303 217.451
            1267 CA ASP
                           131
                                    71.442
                                                            1,00 49.58
                                            82.685 218.322
                                                            1.00 53.33
     MOTA
            1276 CA TRP
                           132
                                    67.946
40
     MOTA
            1292
                  CA ASN
                           133
                                    67.568
                                            85.178 221.086
                                                            1.00 59.61
     MOTA
            1303
                  CA ILE
                           134
                                    68.572
                                            88.447 219.305
                                                            1.00 49.82
     MOTA
            1312 CA PHE
                           135
                                    64.945
                                            89.516 219.045
                                                            1.00 54.55
     MOTA
            1324 CA SER
                            136
                                    64.745
                                            89.689 222.796
                                                            1.00 62.26
45
                                            93.119 222.569
                                                            1.00 54.38
     MOTA
            1332 CA LYS
                            137
                                    66.563
                                            95.734 224.504
                                                            1.00 58.71
     MOTA
            1345 CA ASN
                                    64.785
                            138
     MOTA
                                    63.902 98.337 222.050 1.00 51.21
            1356 CA CYS
                            139
```

55

```
MOTA
             1363 CA ASN
                             140
                                     61.704 100.578 224.157
                                                              1.00 50.25
      MOTA
             1374
                   CA ASN
                             141
                                     64.219 103.338 224.061
                                                              1.00 58.30
      MDTA
5
             1385
                   CA SER
                             142
                                     65.158 102.992 220.352
                                                              1.00 52.18
      ATOM
             1393
                  CA PHE
                             143
                                     61.498 103.434 219.645
                                                              1.00 45.87
      MOTA
             1405 CA ALA
                             144
                                     61.829 106.488 221.831
                                                              1.00 47.49
      MOTA
             1411
                  CA GLU
                             145
                                     64.873 108.022 220.050
                                                              1.00 49.94
10
      MOTA
             1421 CA CYS
                             146
                                     62.489 108.299 217.017
                                                              1.00 53.65
      MOTA
                   CA SER
             1428
                             147
                                     60.005 111.197 217.009
                                                              1.00 57.45
      MDTA
             1436 CA SER
                             148
                                     57.163 110.805 214.395
                                                              1.00 59.53
      MOTA
             1444 CA ALA
                             149
                                     54.284 108.658 214.945
                                                              1.00 59.20
15
      MOTA
             1450 CA GLY
                             150
                                     51.655 108.783 212.248
                                                              1.00 61.95
      ATOM
              1455 CA HIS
                             151
                                     48.764 106.945 210.680
                                                              1.00 69,20
       MOTA
              1468 CA GLU
                             152
                                      46.475 107.377 207.803
                                                              1.00 76.61
       ATOM
              1478 CA ALA
                             153
                                      45.813 106.188 204.264
                                                              1.00 78.85
20
       MOTA
              1492 CA SER
                             404
                                      43.875
                                             81.916 155.536
                                                              1.00 38.68
       MOTA
              1495 CA GLU
                             405
                                      41.939
                                              78.947 156.798
                                                              1.00 43.57
       MOTA
              1505 CA TYR
                             406
                                      44.549
                                              78.204 159.422
                                                              1.00 41.97
25
       MOTA
             1519 CA CYS
                             407
                                      43.251
                                              81.196 161.434
                                                              1.00 36.99
       MOTA
              1526 CA SER
                             408
                                              78.917 163.442
                                      41.011
                                                              1.00 34.02
       MOTA
              1534 CA HIS
                             409
                                      43.988
                                              77.122 165.021
                                                              1.00 40.02
       MOTA
              1547 CA MET
                             410
                                      46.155 79.994 165.895
                                                              1.00 36.48
30
       MDTA
              1556
                   CA ILE
                             411
                                      44.849
                                              80.843 169.353
                                                              1.00 37.88
       MDTA
              1565 CA GLY
                             412
                                      44.508
                                              77.692 171.141
                                                              1.00 39.09
       MOTA
              1570 CA SER
                             413
                                      43.275
                                              77.480 174.756
                                                             1.00 45.21
       MOTA
              1578 CA GLY
                             414
                                      46.813
                                              77.142 176.120
                                                              1.00 41.00
35
       MDTA
              1583 CA HIS
                             415
                                      47.238
                                              80.801 175.208
                                                              1.00 38.01
       MOTA
              1596 CA LEU
                             416
                                      44.181
                                              81.554 177.458
                                                              1.00 39.12
       ATOM
              1605 CA GLN
                             417
                                      45.501
                                              79.575 180.501
                                                              1.00 39.79
       MOTA
              1617 CA SER
40
                              418
                                      48.547
                                              81.846 180.317
                                                               1.00 36.95
       MOTA
              1625 CA LEU
                             419
                                      46.482
                                              85.106 180.212
                                                              1.00 33.42
       MDTA
              1634 CA GLN
                              420
                                                               1.00 38.37
                                      44.564
                                              83.708 183.162
       MOTA
              1646 CA ARG
                              421
                                      47.839
                                              82.875 185.041
                                                               1.00 41.95
       MOTA
              1663 CA LEU
                              422
                                              86.514 184.405
                                                               1.00 33.04
                                      49.151
       MOTA
              1672 CA ILE
                              423
                                      45.854
                                              87.864 185.738
                                                               1.00 25.14
```

```
85.549 188.861
                                                             1.00 32.93
      MOTA
            1681 CA ASP
                            424
                                     45.824
                                             86.494 189.904
      MOTA
            1690
                  CA SER
                            425
                                     49.394
                                                              1.00 30.92
5
                                             90.276 189.999
                                                              1.00 29.47
                            426
      MOTA
            1698
                  CA GLN
                                     48.776
                                     48.379 91.899 193.372
                                                              1.00 31.04
            1710
                 CA MET
                            427
      MOTA
                                            93.418 193.655
                                                              1.00 42.94
            1719
                 CA GLU
                            428
                                     44.978
      MDTA
                                     45.965 96.999 194.423
                                                              1.00 48.52
      MOTA
             1729
                 CA THR
                            429
10
                                     43.656 99.720 193.224
                                                              1.00 48.67
      MOTA
             1738
                  CA SER
                            430
                                     46.604 101.354 191.803
                                                              1.00 52.34
      MOTA
             1746 CA CYS
                            431
                                     45.813 102.789 188.396
                                                              1.00 50.77
      MDTA
             1753
                  CA GLN
                            432
15
      MOTA
             1765 CA ILE
                            433
                                     47.826 103.166 185.221
                                                              1.00 51.17
                  CA THR
                                                              1.00 52.62
      MOTA
             1774
                            434
                                     46.923 105.057 182.065
                                                              1.00 50.53
      MDTA
             1783 CA PHE
                            435
                                     46.641 103.293 178.672
                                    45.167 103.968 175.164
                                                              1.00 50.56
      MOTA
             1795
                   CA GIU
                             436
20
                                     42.643 101.271 174.193
                                                              1.00 46.56
      MDTA
             1805 CA PHE
                             437
                                     39.404 101.233 172.116
                                                              1.00 52.14
      MOTA
             1817 CA VAL
                             438
                                     36.342 103.107 173.476
                                                              1.00 59.69
                             439
      MOTA
             1825
                   CA ASP
                                     33.750 100.386 173.243
                  CA GLN
                                                              1.00 63.82
      MOTA
             1834
                             440
25 -
                                                              1.00 64.79
                                      30.764 102.682 173.229
      MOTA
             1846 CA GLU
                             441
                                                              1.00 59.01
                                      31,810 104,569 170,091
      MDTA
             1856 CA GLN
                             442
                                                              1.00 55.78
                                      32.994 101.521 168.139
      MOTA
             1868
                  CA LEU
                             443
30
                                      30.279 99.137 169.086
                                                              1.00 54.82
      MOTA
             1877
                   CA LYS
                             444
                                      30.862 96.668 166.093
                                                              1.00 51.31
      MOTA
             1890
                   CA ASP
                             445
                                      32.144 93.267 167.135
                                                              1.00 47.34
      MOTA
             1899
                   CA PRO
                             446
                                                              1.00 41.42
                                      34.363 92.592 164.111
      MOTA
             1906
                   CA VAL
                             447
35
                                                              1.00 39.44
                                              96.184 163.845
      ATOM
             1914
                   CA CYS
                             448
                                      35.680
                                                              1.00 40.33
                                      36.273 96.949 167.607
      MDTA
             1921
                   CA TYR
                             449
                                                              1.00 30.87
                                              93.580 167.624
      MOTA
             1935
                   CA LEU
                             450
                                      38.224
                                      39.989 94.664 164.462
                                                              1.00 34.49
40
       MOTA
              1944
                   CA LYS
                             451
                                                              1.00 35.38
                                      40.962 97.918 166.163
       MOTA
             1957
                   CA LYS
                             452
                                                              1.00 28.00
                                      41.503 96.317 169.572
       MDTA
              1970
                   CA ALA
                             453
                                                              1.00 32.56
                                      44.012 93.845 168.164
       MOTA
              1976
                   CA PHE
                             454
45
                                      46.236 96.653 166.861
                                                              1.00 40.19
       MOTA
              1988
                   CA LEU
                             455
                                              98.605 170.114
                                                              1.00 42.45
                                      46.402
       MOTA
              1997
                   CA LEU
                             456
                                      46.966 95.358 172.060
                                                              1.00 38.09
       MOTA
              2006
                   CA VAL
                             457
```

55

```
MOTA
             2014
                   CA GEN
                             458
                                     50.298
                                              95.232 170.161
                                                             1.00 48.02
      MOTA
             2026
                   CA ASP
                             459
                                     51.630
                                              98.529 171.587
                                                              1.00 41.96
5
      MOTA
             2035
                   CA ILE
                             460
                                     50.230
                                              97.686 175.001
                                                              1.00 25.82
      MOTA
             2044
                   CA MET
                             461
                                     52.409
                                              94.542 175.085
                                                              1.00 33.69
      MOTA
             2053
                   CA GUI
                             462
                                     55.702 96.253 173.969
                                                              1.00 45.36
      MOTA
             2063
                   CA ASP
                             463
                                      55.069
                                              99.093 176.249
                                                              1.00 48,41
10
      MOTA
             2072
                   CA THR
                             464
                                     52.771
                                              98.555 179.186
                                                              1.00 41.63
      MOTA
             2081
                   CA MET
                             465
                                      52.743
                                              94.854 180.171
                                                              1.00 42.05
      MOTA
             2090 CA ARG
                             466
                                      56.333
                                              94.928 181.338
                                                              1.00 44.99
      MOTA
             2107 CA PHE
                             467
                                      57.618
                                              91.923 183.501
                                                              1.00 38.72
15
      MOTA
             2119 CA ARG
                             468
                                      61.258 90.894 184.513
                                                              1.00 50.64
      ATOM
             2136 CA ASP
                             469
                                      63.046
                                              88.868 181.895
                                                              1.00 53.46
      MOTA
             2145 CA ASN
                             470
                                      62.896
                                             85.096 181.832
                                                              1.00 52.84
      MDTA
             2156 CA THR
                             471
20
                                      60.211
                                              84.798 184.479
                                                              1.00 47.58
      MOTA
             2165 CA PRO
                             472
                                      57.055
                                              82.652 183.744
                                                              1.00 46.82
      MOTA
             2172 CA ASN
                             473
                                      54.971
                                              85.777 183.172
                                                              1.00 38.75
       MOTA
             2183 CA ALA
                             474
                                      57.431
                                              87.402 180.828
                                                              1.00 37.27
25
       MOTA
             2189 CA ILE
                             475
                                      57.795
                                              84.219 178.807
                                                              1.00 32.85
       MOTA
             2198 CA ALA
                             476
                                      54.021
                                              84.034 178.455
                                                              1.00 27.98
       MOTA
             2204 CA ILE
                             477
                                      53.844
                                              87.522 176.756
                                                              1.00 30.47
       MOTA
              2213 CA VAL
                             478
                                      56.690
                                              86.495 174.201
                                                              1.00 30.40
30
       MOTA
              2221
                   CA GLN
                             479
                                      54.381
                                              83.711 173.325
                                                              1.00 31.18
       MDTA
              2233
                   CA LEU
                             480
                                              86.038 172.893
                                      51.419
                                                              1.00 27.44
       MDTA
              2242 CA GIN
                             481
                                      53.392
                                              88.551 170.821
                                                              1.00 31.59
       MOTA
              2254
                   CA GLU
                             482
                                                             1.00 29.90
                                      54.572
                                              85.683 168.506
35
       MOTA
              2264 CA LEU
                             483
                                      50.900
                                              84.532 168.340
                                                              1.00 22.16
       MOTA
              2273 CA SER
                             484
                                      50.077
                                              88.156 167.451
                                                               1.00 29.37
       MOTA
              2281 CA LEU
                             485
                                      52.410
                                              87.880 164.490
                                                              1.00 32.70
40
       MOTA
              2290 CA ARG
                             486
                                              84.572 163.229
                                                               1.00 30.59
                                      50.865
       MDTA
              2307 CA LEU
                             487
                                      47.396
                                              86.383 163.343
                                                              1.00 31.02
       ATOM
              2316 CA LYS
                             488
                                      48.292
                                              89.399 161.352
                                                               1.00 34.90
       ATOM
              2329 CA SER
                             489
                                      47.728
                                              87.397 158.200
                                                               1.00 35.68
45
       MOTA
              2337 CA CYS
                              490
                                      44.112
                                              87.134 159.417
                                                               1.00 35.03
       MOTA
              2344 CA PHE
                              491
                                      43.520
                                              90.786 159.602
                                                               1.00 42.73
```

50

```
2356 CA THR
      MOTA
                            492
                                     43.145 93.174 156.663
                                                             1.00 56.24
      MDTA
            2365 CA ALA
                            493
                                                             1.00 57.56
                                     44.877 96.467 156.617
      MOTA
            2371 CA ASP
                            494
                                     42.318 99.142 157.302
                                                             1.00 57.16
      MOTA
            2380 CA ALA
                            495
                                     42.063 102.909 156.516
                                                             1.00 62.14
      MOTA
            2386 CA GLU
                            496
                                     45.542 104.198 157.805
                                                              1.00 64.30
      MOTA
            2396 CA GLU
                            497
                                     44.589 107.547 159.382
                                                             1.00 70.67
10
            2406 CA HIS
                            498
      MOTA
                                     41.577 105.660 160.752
                                                              1.00 75.16
            2419 CA ASP
                            499
      MDTA
                                     44.400 103,329 161.968
                                                             1.00 72.11
      MOTA
            2428 CA LYS
                            500
                                     44.773 105.712 164.874
                                                             1.00 65.63
15
      MOTA
            2441
                  CA ALA
                            501
                                     41.040 106.628 165.617
                                                             1.00 55.63
            2447 CA CYS
                            502
      MOTA
                                     38.563 105.427 168.336
                                                             1.00 53.33
      MOTA
             2454 CA VAL
                            503
                                     41.155 105.639 171.100
                                                             1.00 52.46
             2462 CA ARG
      MOTA
                            504
                                     39.866 106.944 174.425
                                                             1.00 51.00
20
      MDTA
             2479 CA THR
                            505
                                     42.588 106.921 177.117
                                                             1.00 50.42
      MDTA:
             2488
                  CA PHE
                            506
                                     41.955 104.935 180.331
                                                             1.00 43.64
      ATOM
                  CA TYR
             2500
                            507
                                     43.306 105.190 183.809
                                                             1.00 41.84
      ATOM
             2514
                  CA GLU
                            508
                                     42.651 101.611 185.218
                                                              1.00 39.35
25
                            509
                                                              1.00 41.67
      MOTA
             2524
                  CA THR
                                     44.156 99.096 187.597
      MOTA
             2533
                  CA PRO
                            510
                                     46.618 96.563 185.969
                                                              1.00 39.70
      MOTA
             2540
                  CA LEU
                            511
                                     44.134 93.600 186.416
                                                              1.00 42.72
             2549
      MOTA
                  CA GLN
                            512
                                     41.460
                                             95.632 184.624
                                                              1.00 37.47
      MOTA
             2561
                  CA LEU
                            513
                                     43.895 96.148 181.765
                                                              1.00 27.02
      MDTA
             2570
                  CA LEU
                             514
                                             92.421 181.627
                                                              1.00 24.55
                                     44.769
      MOTA
             2579 CA GLU
                            515
                                             91.327 181.688
                                                              1.00 31.93
                                     41.076
35
      MDTA
             2589
                   CA LYS
                            516
                                     40.903
                                             93.442 178.515
                                                              1.00 31.96
                   CA VAL
                                                              1.00 26.46
      MDTA
             2602
                             517
                                     43.981
                                             91.949 176.692
      MOTA
                  CA LYS
                                     42.528 88.538 177.852
                                                              1.00 30.83
             2610
                             518
40
      ATOM
             2623 CA ASN
                             519
                                     39.106 89.397 176.432
                                                              1.00 37.62
      MOTA
             2634 CA VAL
                             520
                                     40.221 90.488 172.933
                                                              1.00 35.75
      MOTA
             2642 CA PHE
                                                              1.00 25.45
                             521
                                     42.191 87.109 172.662
             2654 CA ASN
                             522
                                                              1.00 26.73
      MOTA
                                     39.292
                                             85.079 173.953
45
      MOTA
             2665
                   CA GLU
                             523
                                     36.870
                                             86.722 171.437
                                                              1.00 33.04
      MOTA
                      THR
                             524
                                                              1.00 29.70
             2675
                   CA
                                     39.260
                                              86.923 168.369
      MDTA
                                     39.423
                                             83.195 169.011
                                                              1.00 28.41
             2684 CA LYS
                             525
```

27

50

```
MOTA
              2697 CA ASN
                             526
                                      35.624 82.663 169.124
                                                              1.00 34.46
       MOTA
                             527
              2708
                   CA LEU
                                      34.816 84.653 166.021
                                                              1.00 28.09
       MOTA
              2717 CA LEU
                             528
                                      37.617
                                              82.902 164.068
                                                              1.00 25.62
       MOTA
              2726 CA ASP
                             529
                                      35.899
                                              79.662 164.964
                                                              1.00 35.00
       MDTA
              2735 CA LYS
                             530
                                      32.659
                                              81.068 163.537
                                                              1.00 37.54
       MOTA
              2748 CA ASP
                              531
                                      34.250 82.150 160.272
                                                              1.00 37.92
10
       MOTA
              2757 CA TRP
                              532
                                      38.053
                                              82.453 159.573
                                                              1.00 40.02
       MOTA
              2773 CA ASN
                              533
                                      38.395 85.446 157.159
                                                              1.00 42.64
       MOTA
                       ILE
                              534
              2784 CA
                                      35.840 87.437 159.393
                                                              1.00 39.27
       MOTA
              2793 CA PHE
                              535
 15
                                      38.635 89.996 159.881
                                                              1.00 40.49
       MOTA
              2805 CA SER
                              536
                                      38.485 90.720 156.238
                                                              1.00 46.56
       MOTA
              2813 CA LYS
                              537
                                      35.844 93.371 157.102
                                                              1.00 45.61
        MOTA
              2826 CA 'ASN
                              538
                                      35.653 96.936 155.760
                                                              1.00 54.39
 20
        MDTA
              2837 CA CYS
                              539
                                      35.786 98.626 159.004
                                                              1.00 50.55
        MOTA
              2844
                    CA ASN
                              540
                                      36.738 101.888 157.286
                                                              1.00 58.49
        MOTA
              2855 CA ALA
                              541
                                      33.143 102.934 157.648
                                                              1.00 61.90
        MOTA
              2861 CA SER
                              542
                                      32.854 101.855 161.326
                                                              1.00 56.62
 25
        MOTA
              2869 CA PHE
                              543
                                      35.911 103.888 162.235
                                                               1.00 55.25
        MOTA
              2881
                    CA ALA
                              544
                                      33.885 107.102 161.741
                                                               1.00 61.94
        MOTA
              2887
                    CA
                        GU
                              545
                                      30.949 107.003 164.203
                                                               1.00 63.83
        MOTA
               2897
                    CA CYS
                              546
                                      33.917 107.418 166.572
                                                               1.00 61.90
 30
        MOTA
               2904
                        SER
                              547
                                      33.751 110.878 168.150
                    CA
                                                               1.00 66.22
        MDTA
               2912 CA ALA
                              548
                                      37.461 110.300 168.811
                                                               1.00 65.91
        END
```

35

40

45

55

SEQUENCE LISTING

_	(1) GENERAL INFORMATION:
5	(1) APPLICANTS: CHIRON CORPORATION AND THE REGENTS OF THE UNIVERSITY OF CALIFORNIA
	(ii) TITLE OF INVENTION: CRYSTALLIZATION OF M-CSF
10	(iii) NUMBER OF SEQUENCES: 10
15	(iv) CORRESPONDENCE ADDRESS: (A) ADDRESSEE: Harshall, O'Toole, Geratein, Hurray & Borun (B) STREET: 6300 Sears Tower, 233 South Wacker Drive (C) CITY: Chicago (D) STATE: Illinois (E) COUNTRY: USA (F) EIP: 60606-6402
20	 (v) COMPUTER READABLE FORM: (A) MEDIUM TYPE: Floppy disk (B) COMPUTER: IBM PC compatible (C) OPERATING SYSTEM: PC-DOS/MS-DOS (D) SOFTWARE: PatentIn Release #1.0, Version #1.25
25	(vi) current application data: (A) application number: (B) Filing date: (C) classification:
	(VII) PRIOR APPLICATION DATA: (A) APPLICATION NUMBER: US 07/896,512 CIP (B) FILING DATE: 09-JUN-1992
30	(viii) ATTORNEY/AGENT INFORMATION: (A) NAME: Clough, David W. (B) REGISTRATION NUMBER: 36,107 (C) REFERENCE/DOCKET NUMBER: 31436
35	(ix) TELECOMMUNICATION INFORMATION: (A) TELEPHONE: 312/474-6300 (B) TELEFAX: 312/474-0448 (C) TELEX: 25-3856
40	(2) INFORMATION FOR SEQ ID NO:1: (1) SEQUENCE CHARACTERISTICS:
	(A) LENGTH: 56 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear
45	(ii) MOLECULE TYPE: DNA
	(xi) SEQUENCE DESCRIPTION: SEQ ID NO:1:
50	AGGTGTCTCA TAGAAAGTTC GGACGCAGGC CTTGTCATGC TCTTCATAAT CCTTGG 56
	(2) INFORMATION FOR SEQ ID NO:2:

29

5	(i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 58 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
	(ii) MOLECULE TYPE: DNA	
10	(xi) SEQUENCE DESCRIPTION: SEQ ID NO:2:	
~	CAGGAGAAAG CTTATGTCTG AATATTGTAG CGCCATGATT GGGAGTGGAG CCCTGCAG	58
	(2) INFORMATION FOR SEQ ID NO:3:	
15	(i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 58 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	· .·
20	(ii) MOLECULE TYPE: DNA	
	•	
	(xi) SEQUENCE DESCRIPTION: SEQ ID NO:3:	
25	CAGGAGAAAG CTTATGTCTG AATATTGTAG CGCCATGATT GGGAGTGGAC ACCTGCAG	- 58
	(2) INFORMATION FOR SEQ ID NO:4:	
30	(i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 58 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
	(ii) MOLECULE TYPE: DNA	
35	(xi) SEQUENCE DESCRIPTION: SEQ ID NO:4:	
	CAGGAGAAAG CTTATGTCTG AATATTGTAG CCACATGATT GGGAGTGGAG CCCTGCAG	50
	(2) INFORMATION FOR SEQ ID NO:5:	58
10	(i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 72 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single	
5	(ii) MOLECULE TYPE: DNA	
	(x1) SEQUENCE DESCRIPTION: SEQ ID NO:5:	
o	ACCACADAGE TTATCTCTCA ATATTCTACE CACATGATTC GGAGTGGACA CCTCCAGTCT	60
	CTGGCTCGGC TG	
		72

	(2) INFORMATION FOR SEQ ID NO:6:	
.	(i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 85 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
10	(ii) MOLECULE TYPE: DNA	
	(xi) SEQUENCE DESCRIPTION: SEQ ID NO:6:	
	GGACGCAGGC CITGTCATGC TCTTCATAAT CCTTGGTGAA GCAGCTCTTC AGCCTCAAAG	60
15	AGAGTTCCTG CAGCTGTTTA ATGGC 85	
	(2) INFORMATION FOR SEQ ID NO:7:	•
20	(i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 76 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
25	(11) MOLECULE TYPE: DNA	•
	(xi) SEQUENCE DESCRIPTION: SEQ ID NO:7:	
	CATGACAAGG CCTGCGTCCG AACTTTCTAT GAGACACCTC TCCAGTTGCT GGCGAAGGTC	60
30	AAGGCTGTCT TTAATG 76	
	(2) INFORMATION FOR SEQ ID NO:8:	
35	(i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 27 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
40	(ii) MOLECULE TYPE: DNA	
	(xi) SEQUENCE DESCRIPTION: SEQ ID NO:8:	
	GGATCAGGAT CCCTCGGACT GCCTCTC	. 27
45	(2) INFORMATION FOR SEQ ID NO:9:	
50	(i) SEQUENCE CHARACTERISTICS: (A) LENGTH: 26 base pairs (B) TYPE: nucleic acid (C) STRANDEDNESS: single (D) TOPOLOGY: linear	
	(ii) MOLECULE TYPE: DNA	

(xi) SEQUENCE DESCRIPTION: SEQ ID NO:9:

GCGTACCATG GGCCCAGGAG TTCTGC

26

- (2) INFORMATION FOR SEQ ID NO:10:
 - (i) SEQUENCE CHARACTERISTICS:
 - (A) LENGTH: 33 base pairs
 - (B) TYPE: nucleic acid (C) STRANDEDNESS: single
 - (D) TOPOLOGY: linear
 - (ii) MOLECULE TYPE: DNA

15

10

5

(xi) SEQUENCE DESCRIPTION: SEO ID NO: 10:

AGTCGAGGAT CCTCAATCCG GGGGATGCGT GTG

3:3

20

Claims

- A composition comprising crystalline M-CSF, the M-CSF being a dimer of two M-CSF polypeptide monomers, wherein each monomer has substantially the same amino acid sequence as a polypeptide selected from the group consisting of mature M-CSFα polypeptide, mature M-CSFβ polypeptide, and mature M-CSFγ polypeptide.
 - 2. The composition of claim 1 wherein the crystalline dimer has P2₁2₁2₁ symmetry.

30

- 3. The composition of claim 2 wherein the M-CSF is biologically active when resolubilized.
- 4. The composition of claim 1 wherein the two M-CSF polypeptide monomers are disulfide-linked to each other.
- 35 5. The composition of claim 1 wherein each of the two M-CSF polypeptide monomers is an M-CSF α polypeptide.
 - 6. The composition of claim 1 wherein each monomer has an amino acid sequence between about 145 and about 156 amino acids in length, and wherein said sequence has a carboxylic acid terminus located between about 145 and 156 residues measured from the amino terminus of a mature M-CSF polypeptide.

40

- The composition of claim 6 wherein at least one of the polypeptides contains residues 4 to 158 of mature M-CSFα polypeptide.
- 8. The composition of claim 1 wherein at least one of the M-CSF polypeptide monomers is unglycosylated.

15

- A method of crystallizing a M-CSFα dimer having two polypeptide monomers, each containing 146 to 162 amino acids of mature M-CSFα, the method comprising the steps of:
 - a) mixing a solution of the M-CSF α dimer with a precipitant, whereby an M-CSF α mixture is formed;
 - b) precipitating crystalline M-CSF α dimer; and
 - c) isolating the crystalline M-CSFα dimer.

10. The method of claim 9 wherein the M-CSF α dimer is produced from recombinant M-CSF α isolated from a bacterial cell.

55

- 11. The method of claim 9 wherein the precipitant comprises polyethylene glycol.
- 12. The method of claim 9 wherein at least one of the polypeptide monomers is unglycosylated.

- 13. The method of claim 9 wherein at least one of the polypeptide monomers contains residues 4 to 158 of mature M-CSFα polypeptide.
- 14. The method of claim 9 wherein the crystalline M-CSF α dimer is biologically active when resolubilized.
- 15. The method of claim 9 wherein the step of precipitating the M-CSF α crystal comprises equilibrating the M-CSF α dimer mixture with a precipitant solution.
- The method of claim 15 wherein the precipitant solution comprises a higher concentration of precipitant than the
 M-CSFα mixture.
 - 17. The method of claim 15 wherein the step of equilibrating comprises applying the M-CSFα mixture to a surface and equilibrating the applied M-CSFα mixture with a reservoir of the precipitant solution.
- 18. The method of claim 9 wherein the step of precipitating the crystalline M-CSFα dimer comprises a step of applying M-CSFα seed crystals to the M-CSFα dimer mixture.
 - 19. The method of claim 9 wherein the step of precipitating the M-CSFα dimer comprises the step of altering the temperature of the M-CSFα mixture.
 - 20. Crystalline M-CSFα (4-158)

5

20

35

45

50

- 21. The crystalline M-CSF of claim 20 comprising crystals of M-CSF having measurements of at least about 0.3 mm in at least two dimensions.
- 22. The crystalline M-CSF of claim 20 which is capable of diffracting X-ray radiation to produce a diffraction pattern representing the three-dimensional structure of the M-CSF.
- 23. The crystalline M-CSF of claim 22 wherein the diffraction pattern represents the three-dimensional structure of a portion of the M-CSFα (4-158) to an accuracy of at least five angstroms, wherein the portion includes amino acid residues between residue 4 and about residue 153.
 - 24. A method for determining a three-dimensional structure of an M-CSF dimer containing two M-CSF polypeptide monomers, each having between about 146 and 162 amino acids beginning from about residue 1 to about residue 5 of mature M-CSF, the method comprising:
 - a) crystallizing the M-CSF dimer;
 - b) irradiating the crystalline M-CSF to obtain a diffraction pattern characteristic of the crystalline M-CSF; and
 - c) transforming the diffraction pattern into the three-dimensional structure of the M-CSF dimer.
 - 25. A method of using a three-dimensional structure of M-CSF derived from an M-CSF crystal wherein the three-dimensional structure of M-CSF includes an M-CSF receptor-binding region, the method comprising identifying compounds having structures that interact with a receptor-binding region of the three-dimensional structure of M-CSF and function as an M-CSF agonist or antagonist.
 - 26. The method of claim 25 wherein the three-dimensional structure of M-CSF includes alpha-carbon coordinates substantially the same as those of the structural information presented in Appendix 1.
 - 27. A method of identifying M-CSF agonists or antagonists, the method comprising the following steps:
 - a) crystallizing an M-CSF dimer to form at least one M-CSF crystal, the M-CSF dimer containing a group of amino acid residues defining an M-CSF receptor-binding region;
 - b) irradiating an M-CSF crystal produced by the procedure of step (a) to obtain a diffraction pattern of the M-CSF crystal;
 - c) determining the three-dimensional structure of M-CSF from the diffraction pattern, the three-dimensional structure including an M-CSF receptor-binding region; and
 - d) identifying an M-CSF agonist or antagonist compound having a three-dimensional structure that functionally duplicates essential M-CSF receptor-binding solvent accessible residues present in the three-dimensional

structure of the M-CSF receptor-binding region, said M-CSF agonist or antagonist having altered signal transdiction capacity to M-CSF responsive cells.

28. The method of claim 27 wherein said solvent accessible residues do not participate in formation of the M-CSF dimer interface.

Patentansprüche

5

25

35

45

- Zusammensetzung, umfassend kristallines M-CSF, wobei M-CSF ein Dimer aus zwei M-CSF-Polypeptidmonomeren ist, in dem jedes Monomer im Wesentlichen die gleiche Aminosäuresequenz wie ein Polypeptid hat, das aus der Gruppe, bestehend aus reifem M-CSFα-Polypeptid, reifem M-CSFβ-Polypetid und reifem M-CSFγ-Polypeptid, ausgewählt wird.
 - 2. Zusammensetzung nach Anspruch 1, worin das kristalline Dimer P2₁2₁2₁-Symmetrie aufweist.

3. Zusammensetzung nach Anspruch 2, worin M-CSF biologisch aktiv ist, wenn es wieder in Lösung gebracht wird.

- 4. Zusammensetzung nach Anspruch 1, worin die zwei M-CSF-Polypeptidmonomeren disulfidverbrückt sind.
- 20 5. Zusammensetzung nach Anspruch 1, worin jedes der zwei M-CSF-Polypeptidmonomeren ein M-CSFα-Polypeptid ist
 - 6. Zusammensetzung nach Anspruch 1, worin jedes Monomer eine Aminosäuresequenz einer Länge zwischen etwa 145 und etwa 156 Aminosäuren aufweist und wobei die Sequenz einen Carboxytterminus aufweist, der sich etwa zwischen den Resten 145 und 156 befindet, gemessen von dem Aminoterminus des reifen M-CSF-Polypeptids.
 - Zusammensetzung nach Anspruch 6, worin mindestens eines der Polypeptide die Reste 4 bis 158 des reifen M-CSFα-Polypeptide enthält.
- 30 8. Zusammensetzung nach Anspruch 1, wobei mindestens eines der M-CSF-Polypeptidmonomeren nicht glykosyliert ist.
 - 9. Verfahren zum Kristallisieren eines M-CSFα-Dimers mit zwei Polypeptidmonomeren, wobei jedes 146 bis 162 Aminosäuren des reifen M-CSFα enthält und das Verfahren die folgenden Schritte umfasst:
 - a) Mischen einer Lösung des M-CSF α -Dimers mit einem Fällungsmittel, wobei ein M-CSF α -Gemisch gebildet wird.
 - b) Ausfällen des kristallinen M-CSFα-Dimers und
 - c) Isolieren des kristallinen M-CSFα-Dimers.
 - 10. Verfahren nach Anspruch 9, worin das Fällungsmittel Polyethylenglycol umfasst.
 - Verfahren nach Anspruch 9, worin das M-CSFα-Dimer aus rekombinantem, aus einer Bakterienzelle isoliertem M-CSFα hergestellt wird.
 - 12. Verfahren nach Anspruch 9, worin mindestens eines der Polypeptidmonomere nicht glykosyliert ist.
 - Verfahren nach Anspruch 9, worin mindestens eines der Polypeptidmonomere die Reste 4 bis 158 des reifen M-CSFα-Polypeptids enthält.
 - 14. Verfahren nach Anspruch 9, worin das kristalline M-CSFα-Dimer biologisch aktiv ist; wenn es wieder in Lösung gebracht wird.
- 15. Verfahren nach Anspruch 9, worin der Fällungsschritt des M-CSFα-Kristalls das Equilibrieren des M-CSFα-Dimergemisches mit einer Fällungsmittel-Lösung umfasst.
 - 16. Verfahren nach Anspruch 15, in dem die Lösung des Fällungsmittels eine h\u00f6here F\u00e4llungsmittelkonzentration umfasst als das M-CSFα-Gemisch.

- 17. Verfahren nach Anspruch 15, in dem der Equilibrierungsschritt das Auftragen des M-CSFα-Gemisches auf eine Oberfläche und das Equilibrieren des aufgetragenen M-CSFα-Gemisches mit einem Reservoir der Lösung des Fällungsmittels umfasst.
- Verfahren nach Anspruch 9, in dem der Fällungsschritt des kristallinen M-CSFα-Dimers einen Schritt umfasst, bei dem M-CSFα-Impfkristalle dem M-CSFα-Dimergemisch zugeführt werden.
 - Verfahren nach Anspruch 9, in dem der Fällungsschritt des M-CSFα-Dimers einen Schritt umfasst, bei dem die Temperatur des M-CSFα-Gemisches verändert wird.
 - 20. Kristallines M-CSF α (4-158).

10

15

25

30

35 '

- 21. Kristallines M-CSF nach Anspruch 20, umfassend M-CSF-Kristalle mit Abmessungen von mindestens etwa 0,3 mm in mindestens zwei Dimensionen.
- 22. Kristallines M-CSF nach Anspruch 20, das in der Lage ist, Röntgenstrahlung zu beugen, um ein Beugungsbild zu erzeugen, das die dreidimensionale Struktur von M-CSF repräsentiert.
- 23. Kristallines M-CSF nach Anspruch 22, wobei das Beugungsbild die dreidimensionale Struktur eines Teils von M-CSFα (4-158) mit einer Genauigkeit von mindestens 5 Angström repräsentiert, und dieser Teil Aminosäurereste zwischen Rest 4 und etwa Rest 153 umfasst.
 - 24. Verfahren zur Bestimmung einer dreidimensionalen Struktur eines M-CSF-Dimers mit zwei M-CSF-Polypeptidmonomeren, wobei jedes etwa zwischen 146 und 162 Aminosäuren aufweist, ausgehend von etwa Rest 1 bis etwa Rest 5 des reifen M-CSF, das Verfahren umfassend:
 - a) Kristallisieren des M-CSF-Dimers;
 - b) Bestrahlen von kristallinem M-CSF, um ein für kristallines M-CSF charakteristisches Beugungsbild zu erhalten und
 - c) Transformieren des Beugungsbildes in eine dreidimensionale Struktur des M-CSF-Dimers.
 - 25. Verfahren zur Verwendung einer aus einem M-CSF-Kristall abgeleiteten dreidimensionalen Struktur von M-CSF, wobei die dreidimensionale Struktur von M-CSF einen M-CSF-Rezeptorbindungsbereich umfasst und das Verfahren die I-dentifizierung von Verbindungen umfasst, die Strukturen haben, die mit einem Rezeptorbindungsbereich der dreidimensionalen Struktur von M-CSF wechselwirken und als M-CSF-Agonist oder -Antagonist fungieren.
 - 26. Verfahren nach Anspruch 25, wobei die dreidimensionale Struktur von M-CSF alpha-Kohlenstoffkoordinaten umfasst, die im Wesentlichen die gleichen sind wie diejenigen der Strukturinformationen in Anhang 1.
- 40 27. Verfahren zur Identifizierung von M-CSF-Agonisten oder -Antagonisten, wobei das Verfahren die folgenden Schritte umfasst:
 - a) Kristallisieren eines M-CSF-Dimers unter Bildung von mindestens einem M-CSF-Kristall, wobei das M-CSF-Dimer eine Gruppe von Aminosäureresten enthält, die einen M-CSF-Rezeptorbindungsbereich definieren;
 - b) Bestrahlen eines durch das Verfahren von Schritt (a) hergestellten M-CSF-Kristalls, um ein Beugungsbild des M-CSF-Kristalls zu erhalten;
 - c) Bestimmung der dreidimensionalen Struktur von M-CSF aus dem Beugungsbild, wobei die dreidimensionale Struktur einen M-CSF-Rezeptorbindungsbereich umfasst und
 - d) Identifizierung einer M-CSF-Agonisten- oder Antagonistenverbindung mit einer dreidimensionalen Struktur, die in funktioneller Weise wesentliche, dem Lösungsmittel zugängliche M-CSF-Rezeptorbindungsreste, die in der dreidimensionalen Struktur des M-CSF-Rezeptorbindungsbereichs vorliegen, dupliziert, wobei der M-CSF-Agonist oder -Antagonist ein verändertes Signalübertragungsvermögen gegenüber auf M-CSF ansprechenden Zellen aufweist.
- 28. Verfahren nach Anspruch 27, wobei die dem Lösungsmittel zugänglichen Reste nicht an der Bildung der M-CSF-Dimergrenzfläche teilnehmen.

Revendications

20

- Composition comprenant du M-CSF cristallin, le M-CSF étant un dimère de deux monomères polypeptidiques de M-CSF, où chaque monomère a sensiblement la même séquence d'acides aminés qu'un polypeptide choisi dans le groupe consistant en le polypeptide M-CSFα mature, le polypeptide M-CSFβ mature et le polypeptide M-CSFγ mature.
 - 2. Composition selon la revendication 1 où le dimère cristallin a la symétrie P2₁2₁2₁.
- 10 3. Composition selon la revendication 2 où le M-CSF est biologiquement actif quand il est resolubilisé.
 - Composition selon la revendication 1 où les deux monomères polypeptidiques de M-CSF sont liés l'un à l'autre par une liaison disulfure.
- Composition selon la revendication 1 où chacun des deux monomères polypeptidiques de M-CSF est un polypeptide M-CSFα.
 - 6. Composition selon la revendication 1 où chaque monomère a une séquence d'acides aminés d'une longueur située entre environ 145 et environ 156 acides aminés, et où ladite séquence a une extrémité carboxy-terminale située à environ 145 à 156 résidus en mesurant à partir de l'extrémité amino-terminale d'un polypeptide M-CSF mature.
 - Composition selon la revendication 6 où au moins l'un des polypeptides contient les résidus 4 à 158 du polypeptide M-CSFα mature.
- Composition selon la revendication 1 où au moins l'un des monomères polypeptidiques de M-CSF est non glycosylé.
 - Procédé de cristallisation d'un dimère de M-CSFα ayant deux monomères polypeptidiques contenant chacun 146 à 162 acides aminés du M-CSFα mature, le procédé comprenant les étapes de:
 - a) mélange d'une solution du dimère de M-CSF α avec un précipitant, de sorte qu'il se forme un mélange de M-CSF α ;
 - b) précipitation d'un dimère de M-CSFα cristallin; et
 - c) isolement du dimère de M-CSFa cristallin.
 - 10. Procédé selon la revendication 9 où le précipitant comprend du polyéthylèneglycol.
 - Procédé selon la revendication 9 où le dimère de M-CSFα est produit à partir de M-CSFα recombiné isolé à partir d'une cellule bactérienne.
 - 12. Procédé selon la revendication 9 où au moins l'un des monomères polypeptidiques est non glycosylé.
 - 13. Procédé selon la revendication 9 où au moins l'un des monomères polypeptidiques contient les résidus 4 à 158 du polypeptide M-CSF α mature.
 - 14. Procédé selon la revendication 9 où le dimère de M-CSFα cristallin est biologiquement actif quand il est resolubi-
- 15. Procédé selon la revendication 9 où l'étape de précipitation du cristal de M-CSFα comprend la mise en équilibre du mélange de dimère de M-CSFα avec une solution de précipitant.
 - 16. Procédé selon la revendication 15 où la solution de précipitant comprend une concentration de précipitant plus élevée que le mélange de M-CSFα.
- 755 17. Procédé selon la revendication 15 où l'étape de mise en équilibre comprend l'application du mélange de M-CSFα à une surface et la mise en équilibre du mélange de M-CSFα appliqué avec un réservoir de la solution de précipitant

EP 0 668 914 B1

- 18. Procédé selon la revendication 9 où l'étape de précipitation du dimère de M-CSFα cristallin comprend une étape d'application de germes cristallins de M-CSFα au mélange de dimère de M-CSFα.
- Procédé selon la revendication 9 où l'étape de précipitation du dimère de M-CSFα comprend l'étape de modification de la température du mélange de M-CSFα.
 - 20. M-CSFa (4-158) cristallin.

5

25

30

35

40

- 21. M-CSF cristallin selon la revendication 20 comprenant des cristaux de M-CSF ayant des mesures d'au moins environ 0,3 mm dans au moins deux dimensions.
 - 22. M-CSF cristallin selon la revendication 20 qui est capable de diffracter les rayons X pour produire un diagramme de diffraction représentant la structure tridimensionnelle du M-CSF.
- 23. M-CSF cristallin selon la revendication 22 où le diagramme de diffraction représente la structure tridimensionnelle d'une partie du M-CSFα (4-158) avec une précision d'au moins cinq angströms, où la partie inclut les résidus d'acides aminés entre le résidu 4 et le résidu 153 environ.
 - 24. Procédé de détermination d'une structure tridimensionnelle d'un dimère de M-CSF contenant deux monomères polypeptidiques de M-CSF ayant chacun entre environ 146 et 162 acides aminés depuis les résidus 1 environ à 5 environ du M-CSF mature, le procédé comprenant:
 - a) la cristallisation du dimère de M-CSF;
 - b) l'irradiation du M-CSF cristallin pour obtenir un diagramme de diffraction caractéristique du M-CSF cristallin; et
 - c) la transformation du diagramme de diffraction en la structure tridimensionnelle du dimère de M-CSF.
 - 25. Procédé d'utilisation d'une structure tridimensionnelle de M-CSF issue d'un cristal de M-CSF où la structure tridimensionnelle du M-CSF inclut une région de liaison de récepteur du M-CSF, le procédé comprenant l'identification de composés ayant des structures qui interagissent avec une région de liaison de récepteur de la structure tridimensionnelle de M-CSF et jouent le rôle d'agoniste ou d'antagoniste du M-CSF.
 - 26. Procédé selon la revendication 25 où la structure tridimensionnelle du M-CSF inclut des coordonnées des carbones alpha sensiblement identiques à celles de l'information structurale présentée dans l'annexe 1.
 - 27. Procédé d'identification d'agonistes ou d'antagonistes du M-CSF, le procédé comprenant les étapes suivantes:
 - a) cristallisation d'un dimère de M-CSF pour former au moins un cristal de M-CSF, le dimère de M-CSF contenant un groupe de résidus d'acides aminés définissant une région de liaison de récepteur du M-CSF;
 - b) irradiation d'un cristal de M-CSF produit par le protocole de l'étape (a) pour obtenir un diagramme de diffraction du cristal de M-CSF;
 - c) détermination de la structure tridimensionnelle du M-CSF d'après le diagramme de diffraction, la structure tridimensionnelle incluant une région de liaison de récepteur du M-CSF; et
 - d) identification d'un composé agoniste ou antagoniste du M-CSF ayant une structure tridimensionnelle qui reproduit fonctionnellement les résidus accessibles aux solvants de liaison de récepteur du M-CSF essentiels présents dans la structure tridimensionnelle de la région de liaison de récepteur du M-CSF, ledit agoniste ou antagoniste du M-CSF ayant une capacité de transduction des signaux modifiée à l'égard des cellules sensibles au M-CSF.
- 28. Procédé selon la revendication 27 où lesdits résidus accessibles aux solvants ne participent pas à la formation de l'interface du dimère de M-CSF.

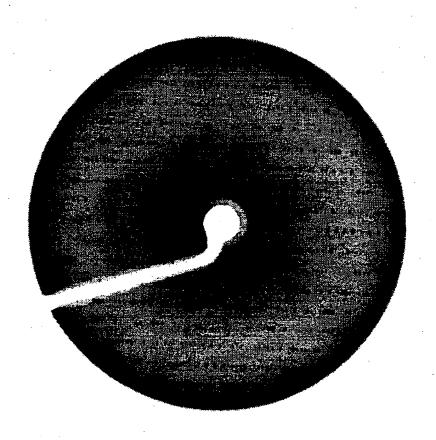


FIG. 1


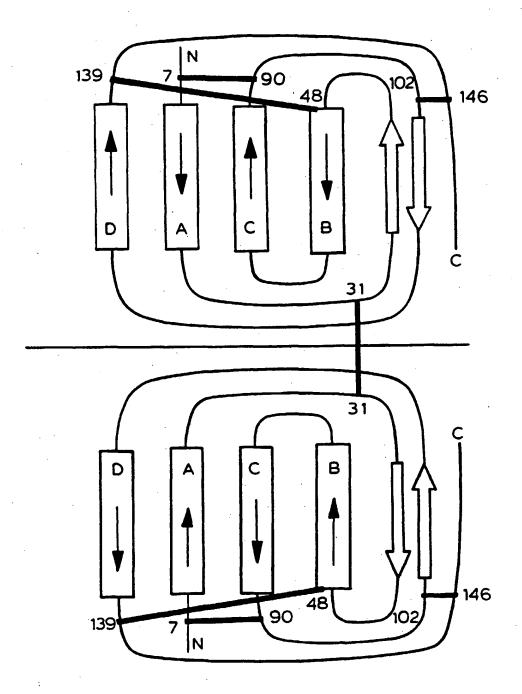


FIG. 2

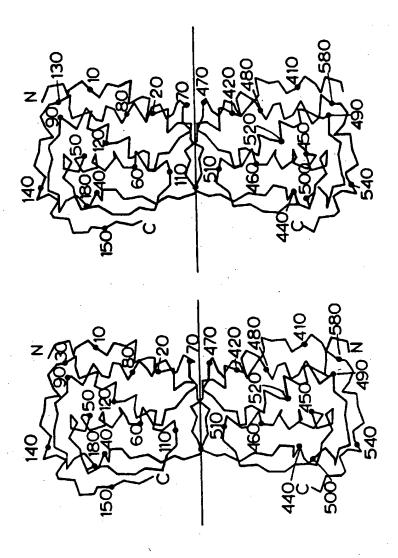

F16.3

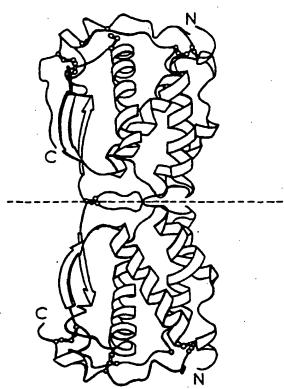


FIG. 4A

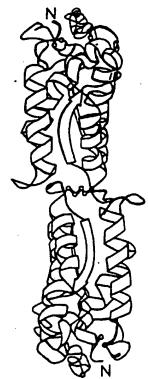
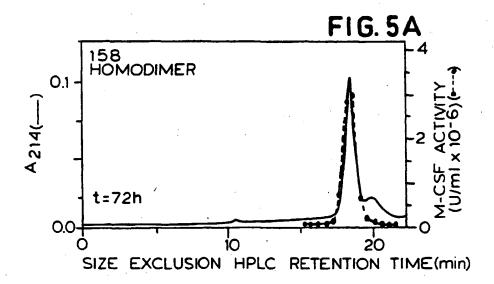
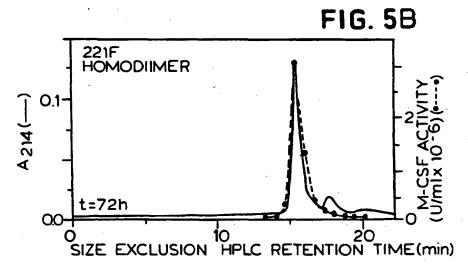
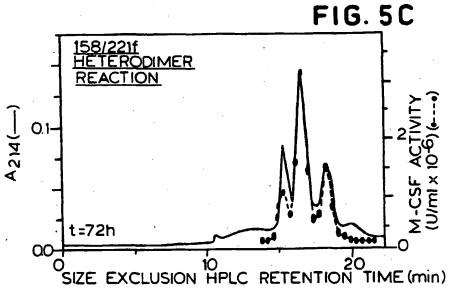


FIG. 48







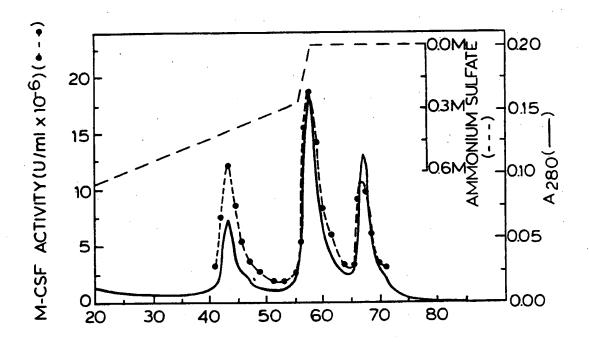


FIG. 6A

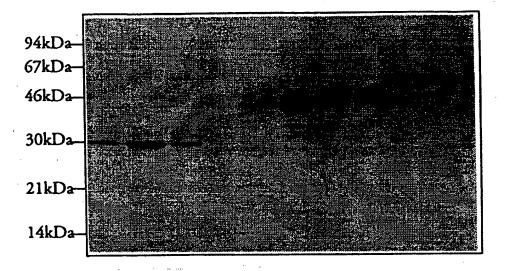


FIG. 6B


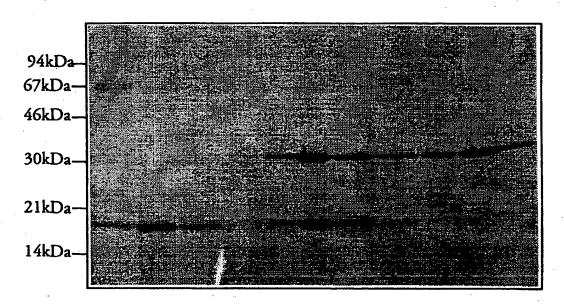


FIG. 6C

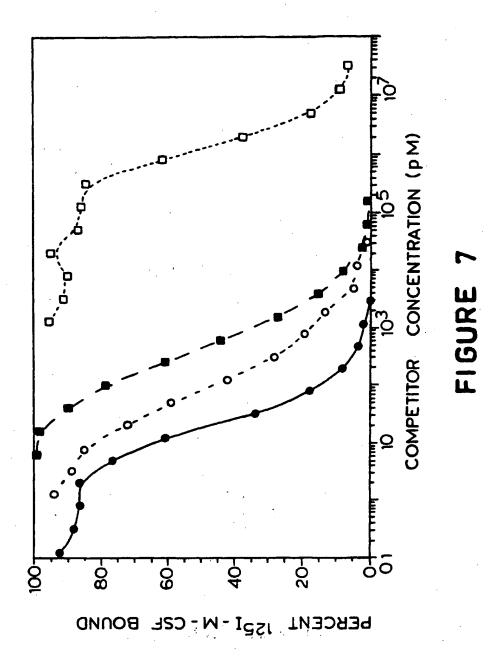

THIS PAGE BLANK (USPTO)

This Page is Inserted by IFW Indexing and Scanning Operations and is not part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:
☐ BLACK BORDERS
☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
☐ FADED TEXT OR DRAWING
□ BLURRED OR ILLEGIBLE TEXT OR DRAWING
☐ SKEWED/SLANTED IMAGES
COLOR OR BLACK AND WHITE PHOTOGRAPHS
GRAY SCALE DOCUMENTS
LINES OR MARKS ON ORIGINAL DOCUMENT
☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
Потигр

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.